


FEATURE GUIDE

WIDEX SOUNDRELAX™

Widex **SoundRelax** in WIDEX **MOMENT™** is the next step for Widex' unique fractal approach, expanding our relaxing sound universe, Zen, with new, softer tones and modulated nature sounds, tailored to the individual user's preferences.

SoundRelax is intended for users for well-being and focus purposes, with or without tinnitus.

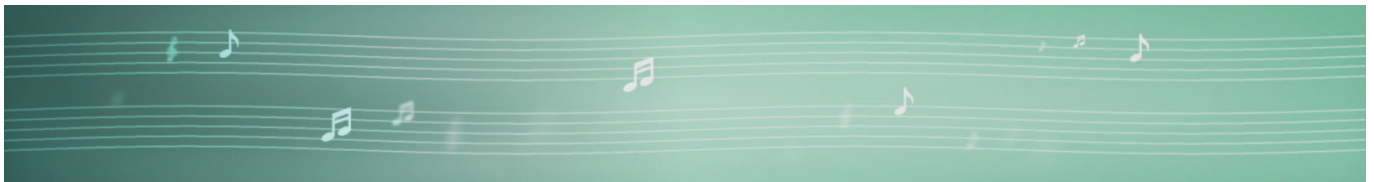

WHAT DOES WIDEX SOUNDRELAX DO?


The new **SoundRelax** sounds include relaxing tones and nature sounds and are designed to both relax and provide acoustic stimulation. **SoundRelax** uses the *principles of fractal mathematics* to create relaxing sounds. Fractal technology ensures that the sounds are being predictable and familiar yet never repeat themselves which in turn helps the brain to de-focus from stress and/or tinnitus.

SoundRelax is built in in the hearing aids to produce the fractal tones. This means that the user does not need to take the hearing aids off to listen to relaxing music and can continue to use the hearing aids and listen to **SoundRelax** while at the same time enjoying the natural sound quality and amplification of WIDEX **MOMENT** hearing aids.

HOW IS WIDEX SOUNDRELAX USED?

SoundRelax is chosen in the software COMPASS GPS 4.4 as either a stand-alone Relax program in the Program Manager or as a Zen+/Relax+ program, which means that the HCP gets three additional program slots to populate within the universe of Zen+/Relax+.


With this update Widex also offers the HCP the option to combine both Zen and **SoundRelax** with Zero-delay PureSound™ as a *program base* option. This is beneficial for the many tinnitus clients that have mild hearing losses, with PureSound, users get the best sound quality which may make them more likely to wear a hearing aid. The program base will as default be the program chosen as P1 in Compass set-up:

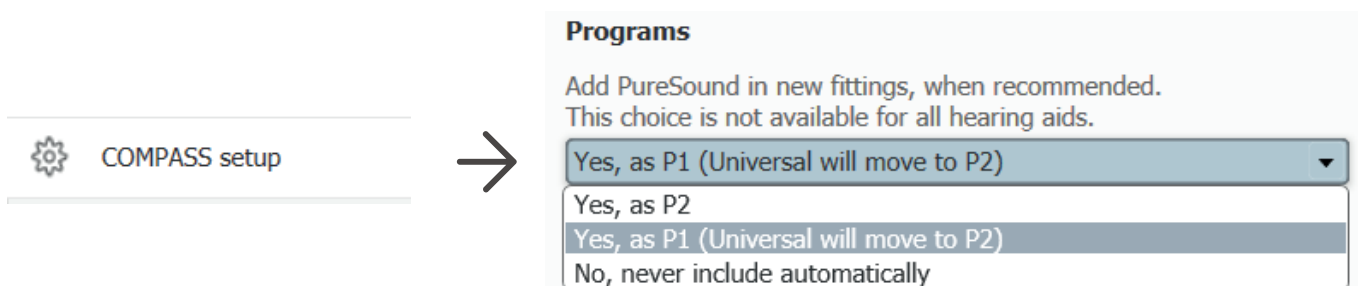


FEATURE GUIDE

WIDEX SOUNDRELAX™

The eight new soundscapes in Widex **SoundRelax** each offers a different soundscape set-up. They are:

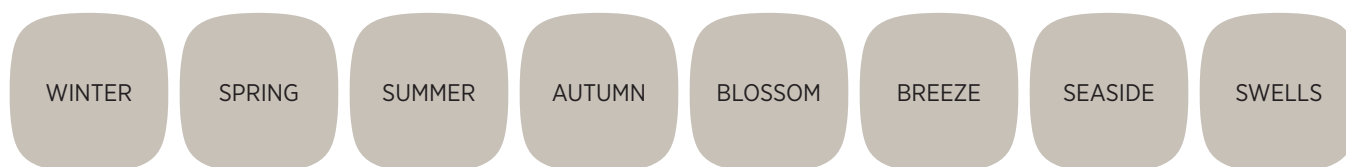

Different settings are available to fine-tune, so that the HCP can find the exact right solution for the user. Choose from:

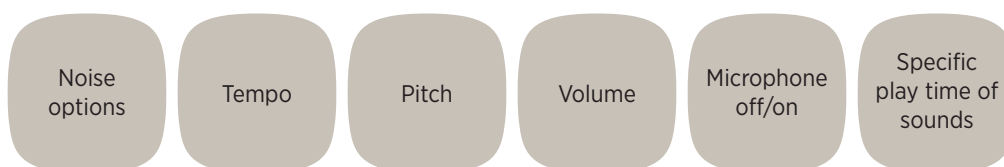

The table below shows an overview of the set-up for each **SoundRelax** style*:

	DEFAULT PITCH			DYNAMIC RANGE		DEFAULT TEMPO			TIMBRE				DEFAULT WAVES	
FRACTAL STYLE	LOW	MEDIUM	HIGH	RESTRICTED	BROAD	SLOW	MEDIUM	FAST	CLEAR	SOFT/PURE	FLUTES	ORGANIC	ON	OFF
SUMMER	■			■			■					■	■	
BREEZE	■			■		■				■				■
SPRING			■	■			■			■			■	
WINTER		■			■		■		■					■
BLOSSOM		■			■		■				■		■	
AUTUMN		■		■			■					■		

*The noise soundscapes do not include tones and are therefore not represented in the table

Users get easy access to the new Relax and Relax+ programs via the program button on the hearing aid, the **MOMENT** app or **TONELINK** app or by using an RC-DEX remote control.

WHAT ARE THE BENEFITS OF WIDEX SOUNDRELAX?

End-user benefits include:

- Ease of use - no interruptions in life
- Never missing out of amplification, feeling safe having all-day use of hearing aids
- Increased use of hearing aids because of the natural sound provided by PureSound

HCP benefits include:

- Offer the best-in-industry solution with unique fractal tones
- Cater to more user's needs with more options for relax and tinnitus handling
- Offer a complete solution with all day use of hearing aids