



UNLOCK THE POTENTIAL: WIDEX SOUND ASSIST™

PRESENTED BY: SUSAN DE BONDT, AU.D.



ABOUT YOUR SPEAKER



Susan de Bondt, AuD
Sr. Clinical Education Specialist
Susan.deBondt@widexsound.com
941-744-6400

SPEAKER DISCLOSURE

SUSAN DEBONDT, AU.D.

Relevant Financial Relationships:

- Presenter is an employee of Widex USA and receives financial compensation.

Relevant Nonfinancial Relationships:

- There are no nonfinancial relationships.

LEARNING OUTCOMES

After this course participants will be able to:

1. List the six functions available in the Widex Sound Assist™ accessory
2. Identify the client needs met by each function of the Widex Sound Assist™
3. Describe the steps of pairing the Widex Sound Assist™ with MOMENT hearing aids.

AGENDA

- Introduction to Widex Sound Assist™
- Review details of each mode of functionality
- Steps for pairing to MOMENT hearing aids
- Steps for pairing to Bluetooth devices
- Quick summary of Sound Assist™ button function
- Access to Sound Assist™ in the MOMENT App

AN ASSIST FOR ANY SITUATION



PARTNER MIC



BLUETOOTH STREAMING



T-COIL



TABLE MIC



HANDS-FREE MIC



REMOTE CONTROL

WHEN TO FIT




- With any Moment™ 2.4 GHz enabled hearing aid
- Provide extra support for:
 - Speech-in-noise complaints
 - Complex listening situations
- Any patient can benefit!



SOUND ASSIST—WHAT'S WHAT?



WHAT'S IN THE BOX?

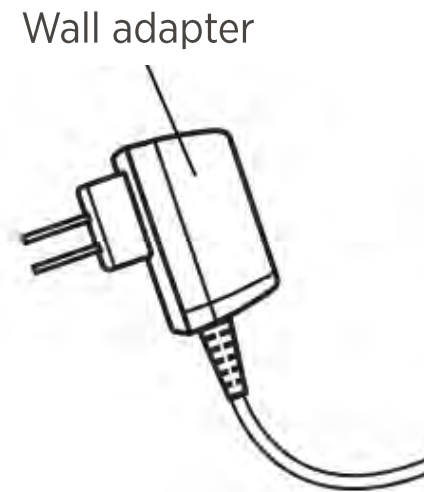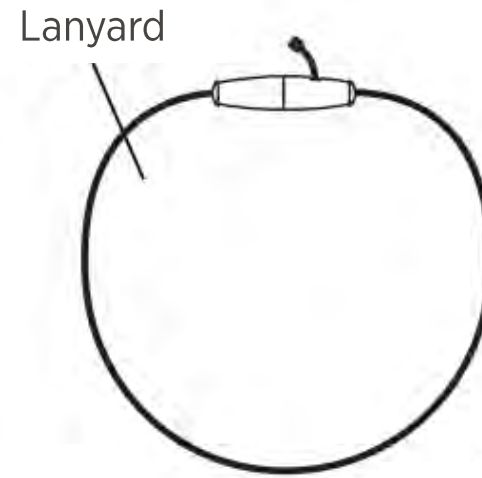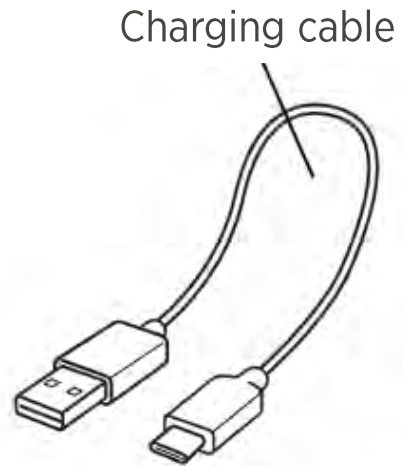

User Manual



SOUND ASSIST STATS

Full charge in

3

hours



Up to*

10

Hours of use



Pair with

1

Set of HAs



Pair with

8

Bluetooth
devices



*Use time will vary based on type of use

AN ASSIST FOR ANY SITUATION



REMOTE MICROPHONE MODE



BLUETOOTH MODE



REMOTE CONTROL



TELECOIL MODE



REMOTE MICROPHONE MODE



PARTNER MIC



TABLE MIC

PARTNER MIC



PARTNER MIC FUNCTION

- 1** Turn Sound Assist on
- 2** Have communication partner wear Sound Assist on lanyard or clip to clothing
- 3** The Sound Assist will automatically go into directional mode facing the speaker up

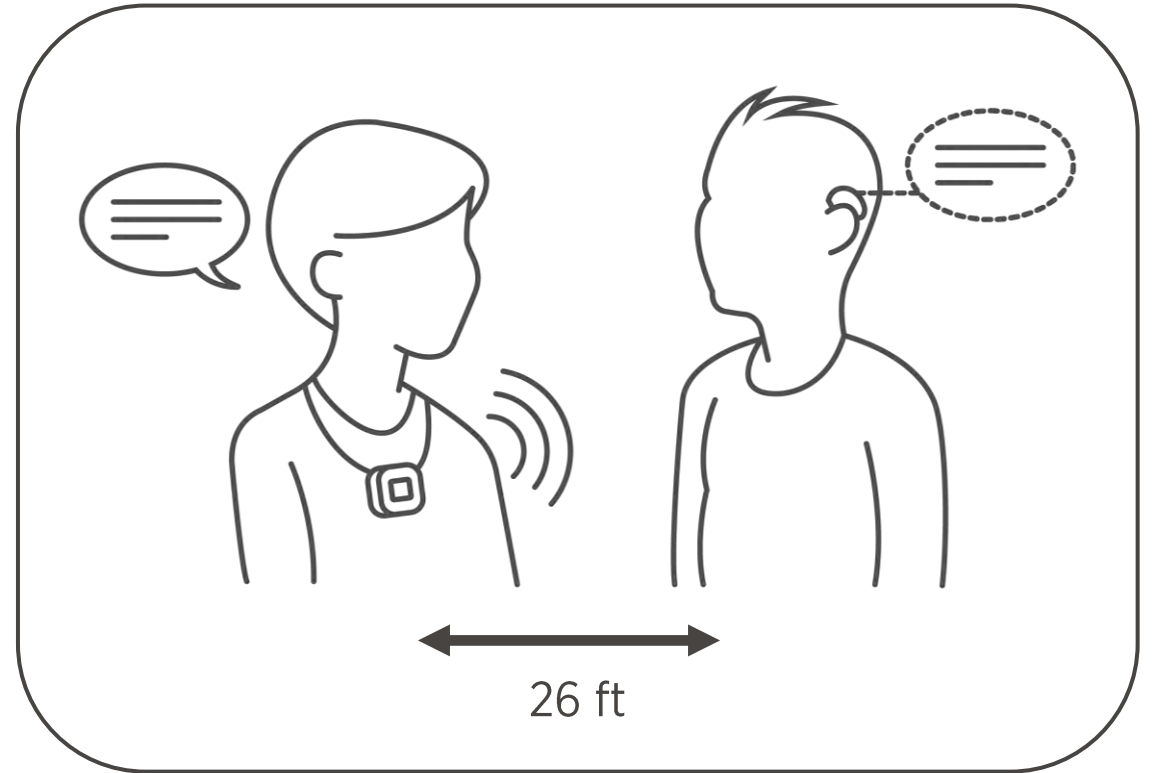


TABLE MIC



TABLE MIC FUNCTION


- 1** Turn Sound Assist on
- 2** Set Sound Assist on table
- 2** *Optional* long press center push-button to lock into directional mode and point in the direction you want it to focus



All mics active
(default)



Listening direction
locked

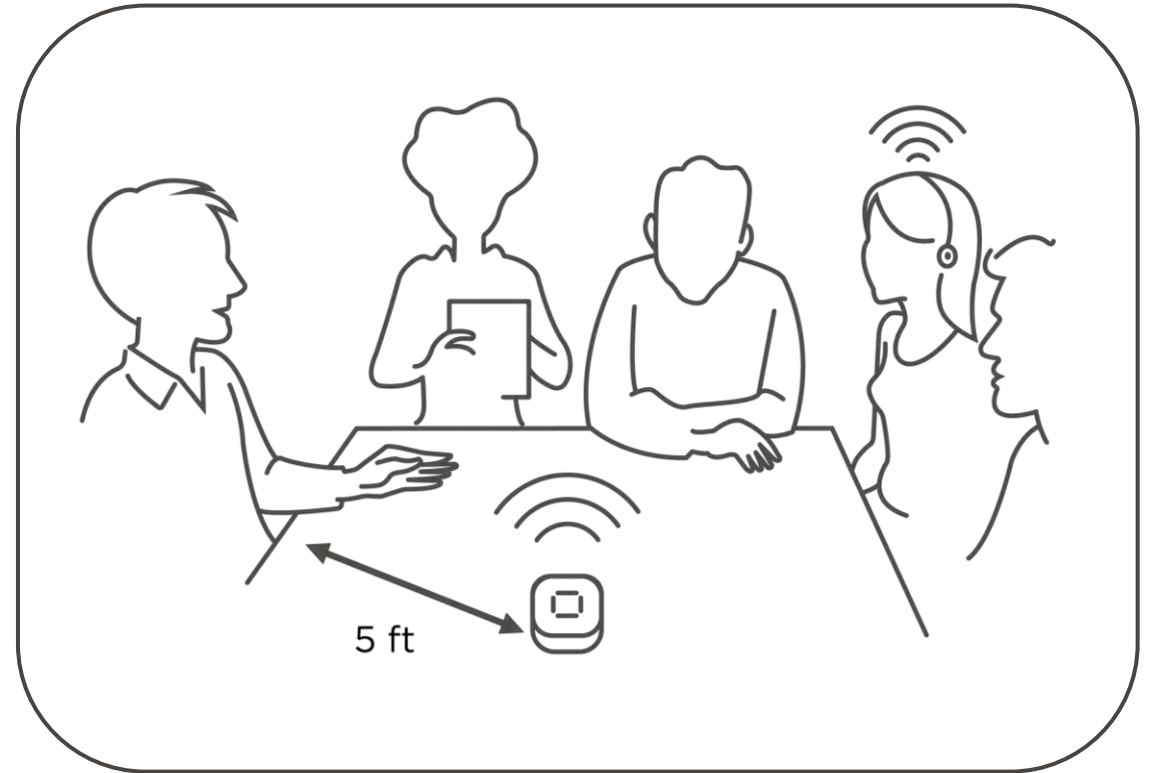

REMOTE MICROPHONE MODE



BLUETOOTH MODE



REMOTE CONTROL



TELECOIL MODE



BLUETOOTH MODE



BLUETOOTH STREAMING



HANDS-FREE MIC


REMOTE CONTROL

BLUETOOTH STREAMING



BLUETOOTH STREAMING FUNCTION

- 1** Pair to any Bluetooth device
- 2** Turn Sound Assist on
- 3** Short press Power button once to toggle to Bluetooth mode
- 4** Start streaming

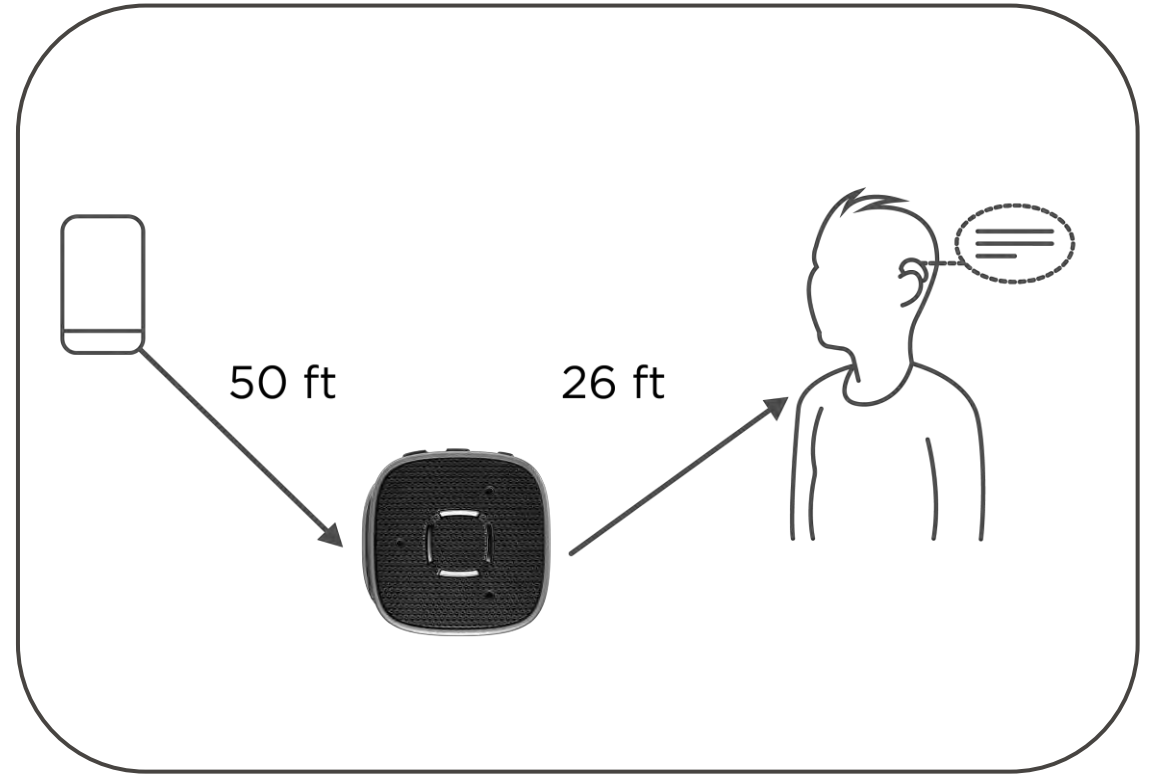


HANDS-FREE MIC



HANDS-FREE MIC FUNCTION


- 1** Pair Sound Assist with any Bluetooth phone
- 2** Wear Sound Assist on lanyard or clip on clothing
- 3** Press center button to answer/end calls
- 4** Long press center button to activate Siri/Voice Assist

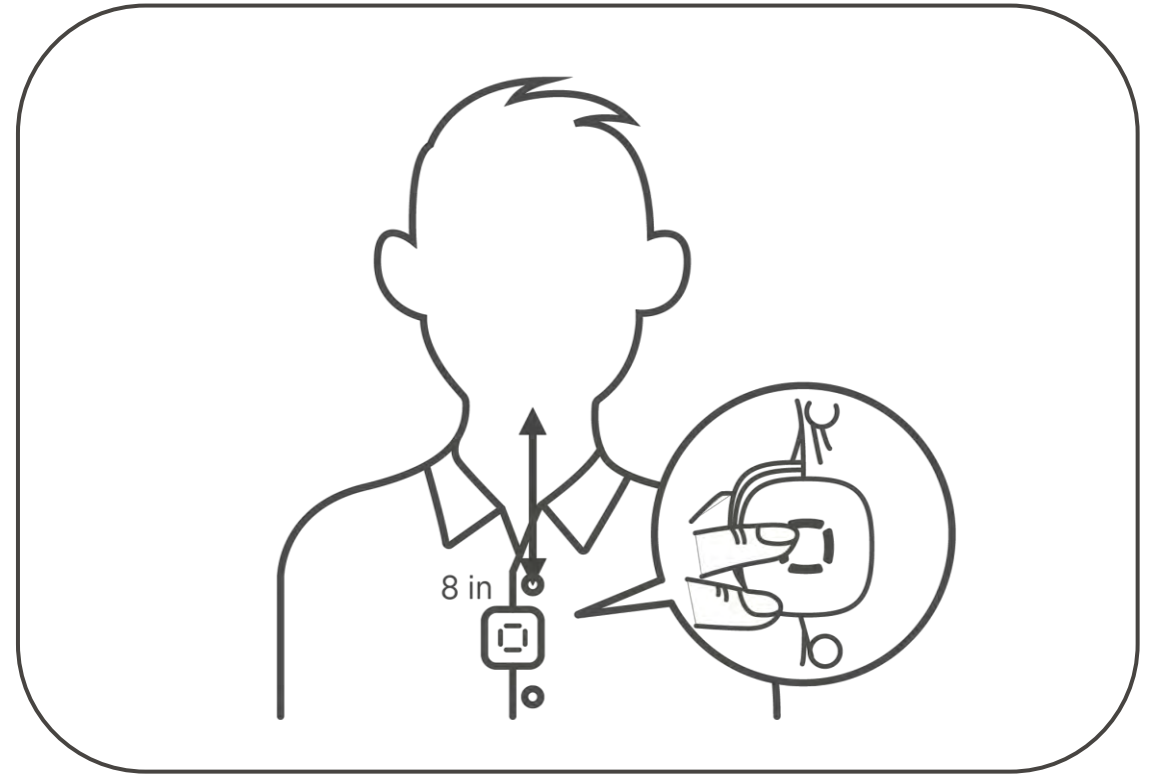

REMOTE CONTROL FUNCTION

Reminders:

- Must be in Bluetooth mode
- Will control HAs when not actively streaming
- If (+) or (-) is pushed while streaming, it will change the volume of the streamed signal

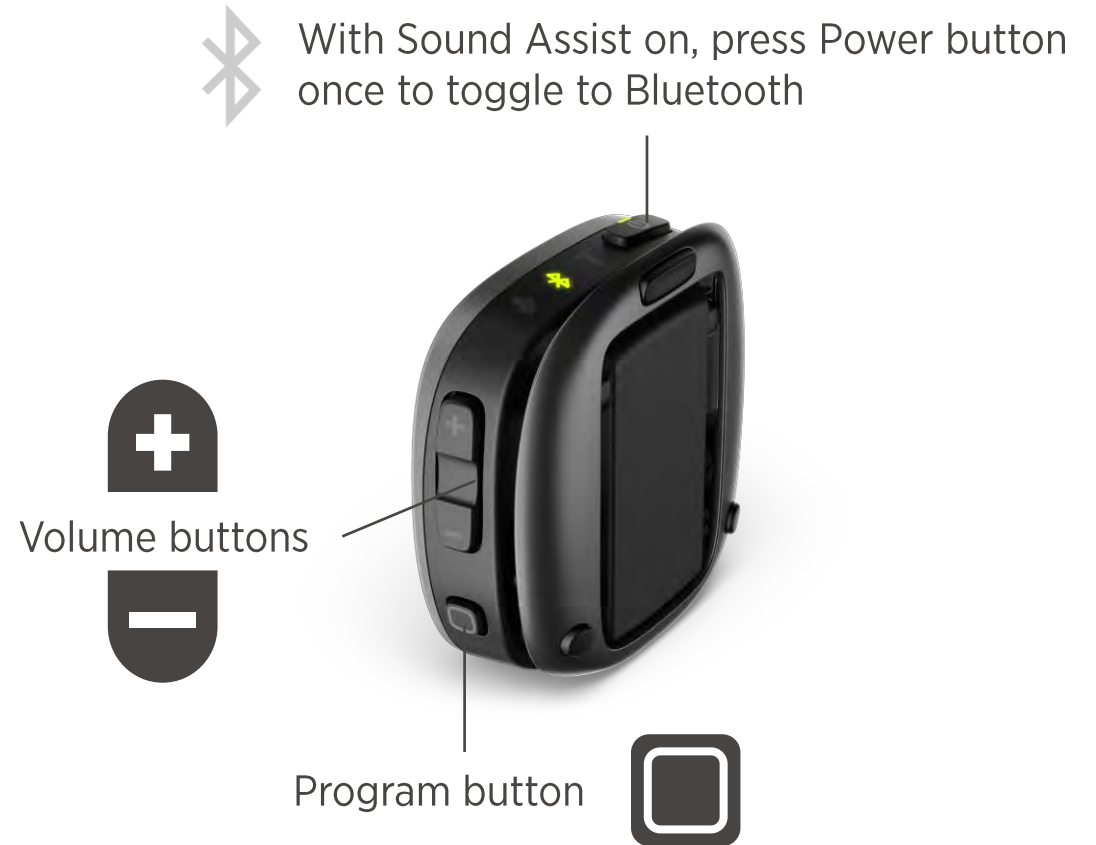

REMOTE MICROPHONE MODE



BLUETOOTH MODE



REMOTE CONTROL



TELECOIL MODE



TELECOIL MODE



T-COIL

T-COIL





T-COIL FUNCTION

- 1** Turn Sound Assist on
- 2** Short press power button twice
- 3** Positioning is important.
Lanyard recommended.
- 4** Use within any telecoil loop system

Power button and indicator light
must be facing up



LANYARD USE RECOMMENDED IN T-COIL MODE

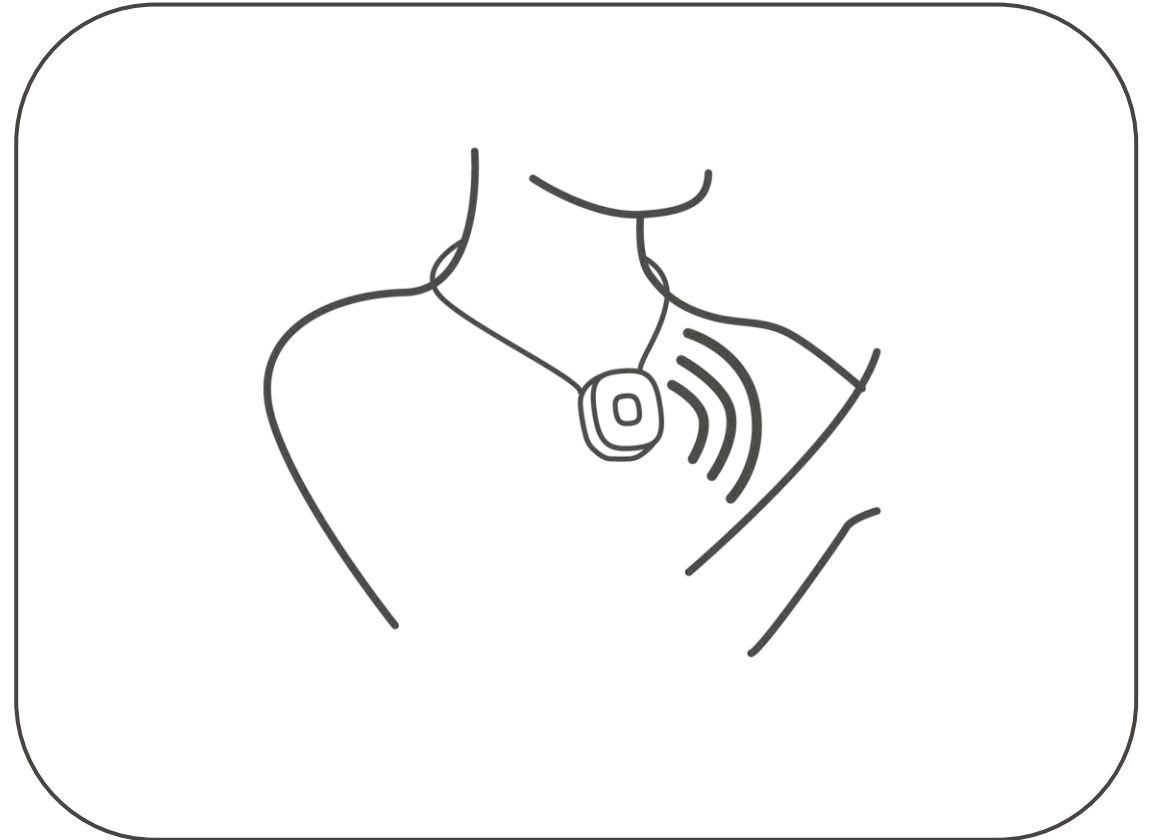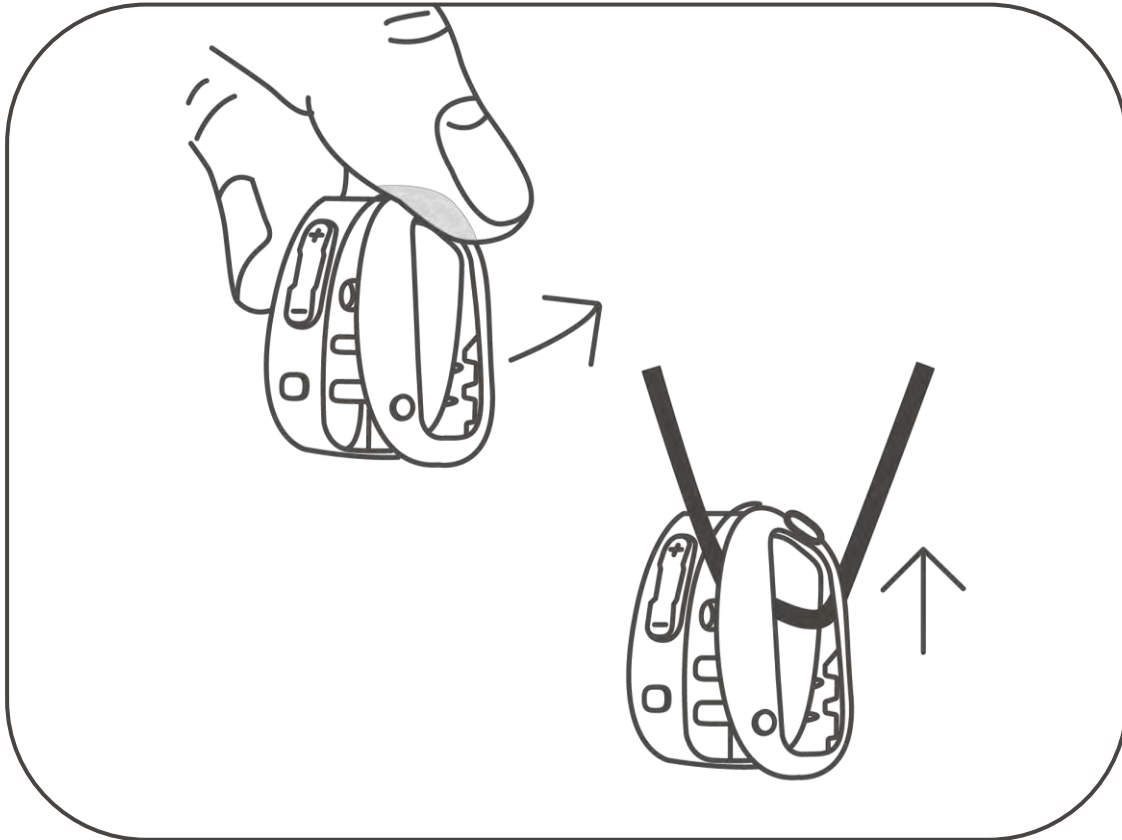



Widex Sound Assist

HEARING AID PAIRING

Press and hold simultaneously
the 'on/off' and '+' buttons for ~5 sec



BASIC FUNCTIONS



Power On/Off

- LONG press (2 sec)

Change Mode

- SHORT press to toggle through modes




Put in Pairing Mode

- With device off
LONG press (3 sec)



Power/Mode button

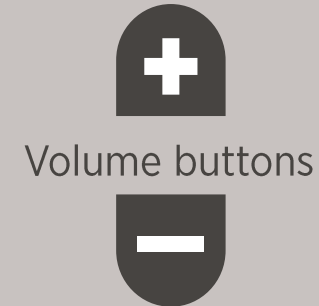

Volume buttons



POWER/MODE BUTTON OVERVIEW



Main Uses:

1. Turns device ON/OFF
2. Toggles between the 3 modes

It is also used for putting the Sound Assist in pairing mode and removing all pairings. With Sound Assist off:



Put in pairing mode



Remove all pairings

CENTER PUSH-BUTTON OVERVIEW



Main Uses:

Pause/resume streaming (regardless of what mode it is in)

- Pause/resume remote mic streaming
- Pause/resume BT streaming
- Pause/resume T coil streaming
- Taking/ending a call

It is also used for activating voice assist when paired to a BT device that has voice assist active (long press)

VOLUME BUTTONS OVERVIEW



Main Uses:

Control volume of incoming signal (regardless of mode)

- VC of remote mic streaming
- VC of BT media streaming
- VC of call streaming
- VC of telecoil streaming
- VC of HA mics (when in BT mode and not streaming)



Long press mutes HA mics (room off feature if pressed while streaming)



Long press unmutes HA mics


PROGRAM BUTTON OVERVIEW



Only Uses:


1. Toggle programs while using as HA remote control
 - Short press for regular programs
 - Long press for Zen/SR+ programs
2. Lock SA buttons (press program button and (+) button simultaneously)

NO PAIRING NEEDED IN COMPASS™ GPS

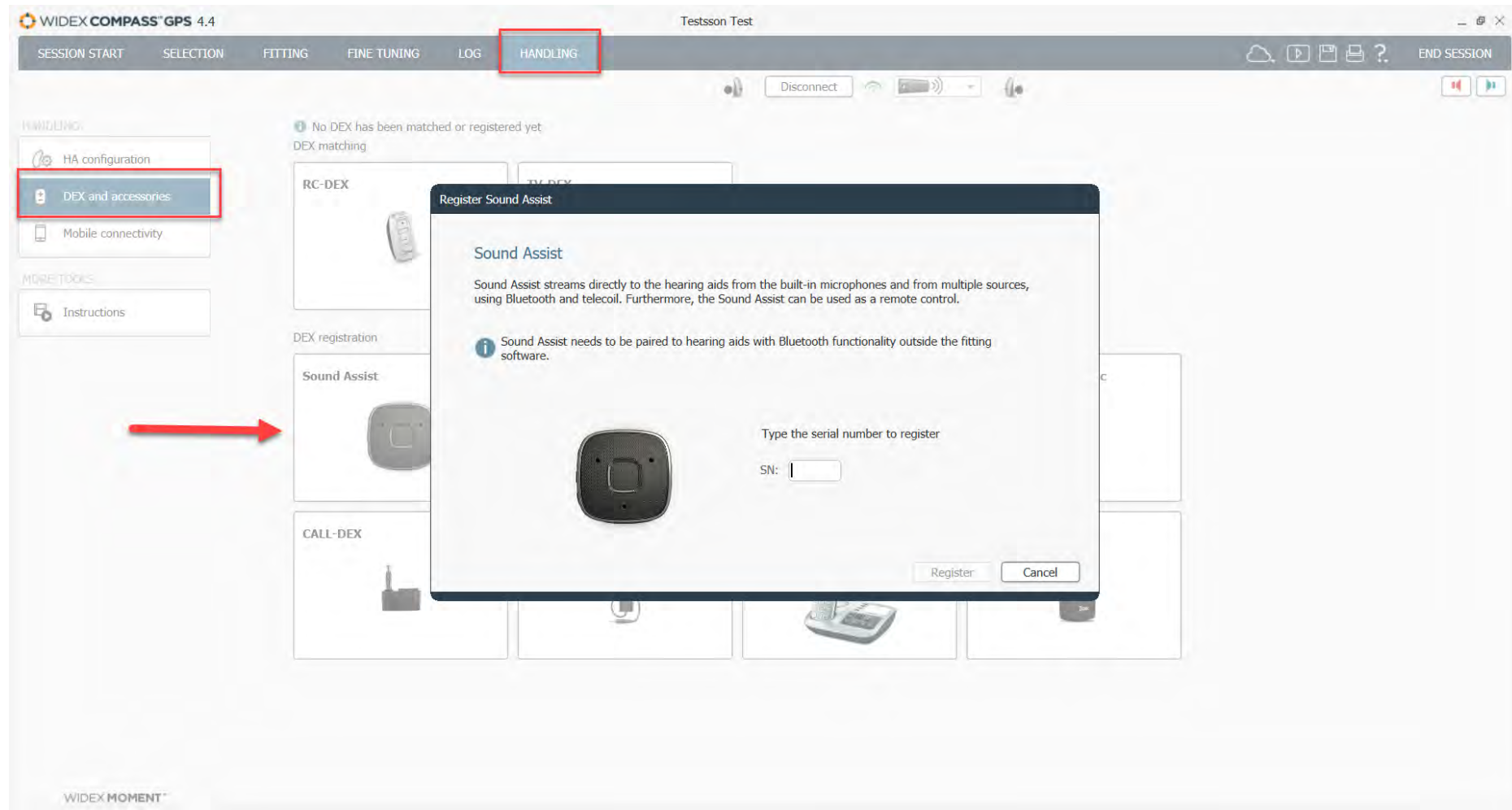


MOMENT™ APP V1.6

- Shows current function
- Sound adjustments available
 - Mute HA mics (room off)
 - Overall streaming volume
 - Sound Mixer
 - L/R Volume
 - Equalizer

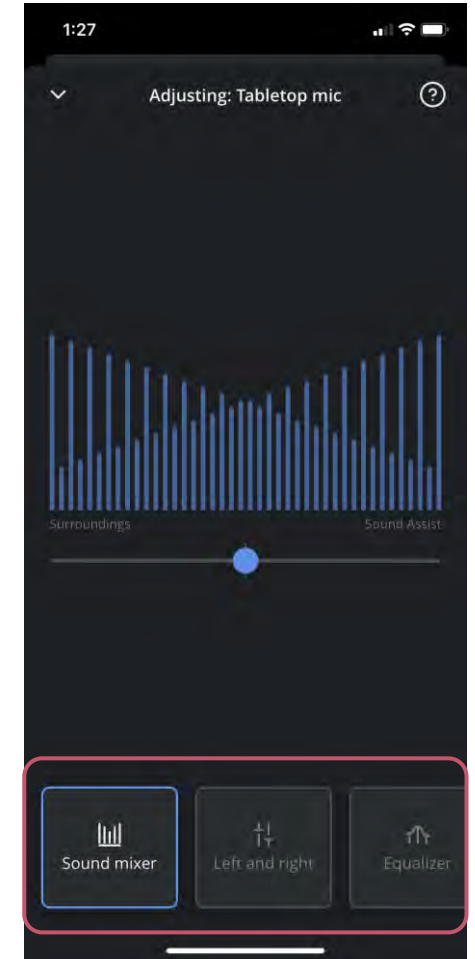

MY GUIDE—SOUND ASSIST

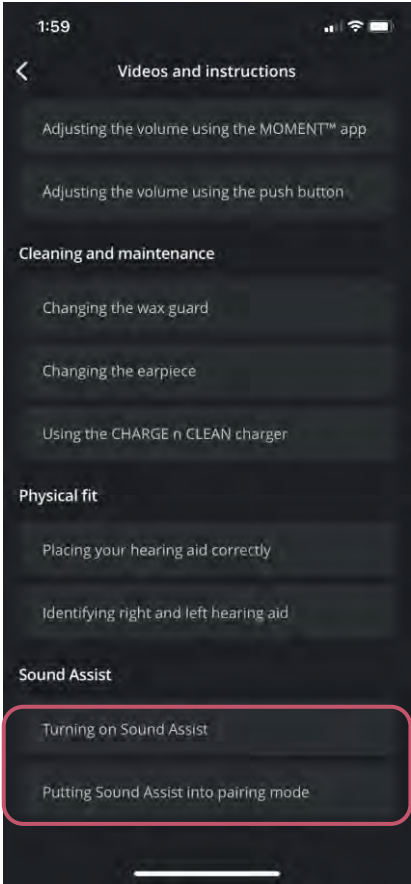

EVERYONE CAN USE AN ASSIST SOMETIMES

- Extremely versatile
- Perfect HA companion
- Excellent price point





SOUND LIKE NO OTHER

Sound Assist video tutorials

Widex.com/en-us

WidexPro.com/en-us





SOUND LIKE NO OTHER

Send your questions to
AQInquiry@widexsound.com