


COMPASS™ GPS FITTING QUICK GUIDE

WIDEX HEARING AIDS

Connect and Initial Setup


1. Choose your programmer and then select **Connect** to detect your Widex hearing aids.

NOTE: If you do not find your preferred programmer in the drop down list, go to Session Start > Compass setup > Interface to configure your programmer options.

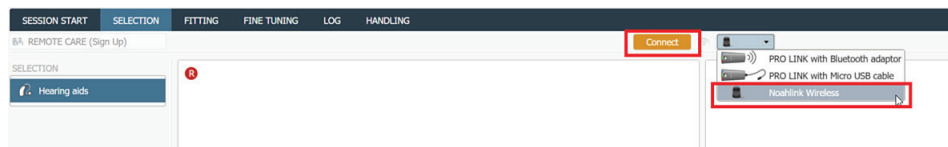


2. Once connected, identify right and left ear by confirming the serial numbers or presenting the musical tone. If applicable, press the button on the device to aid in identification.

TIP: If the devices have been previously connected, confirm what you are seeing on the screen is marked correctly by **selecting the button for right and left**.

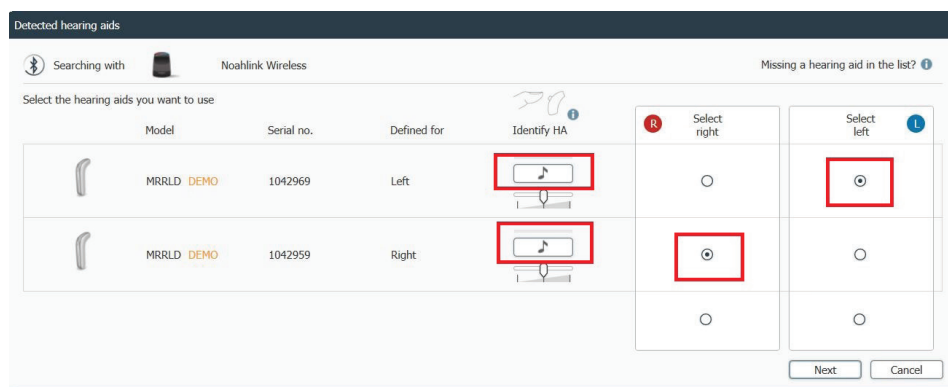

3. Once the aids are connected, the recommended receiver and ear-tip options will be filled in. These selections can be changed through the drop-down menus. The receivers and ear-tips selected should match what is on the aids.

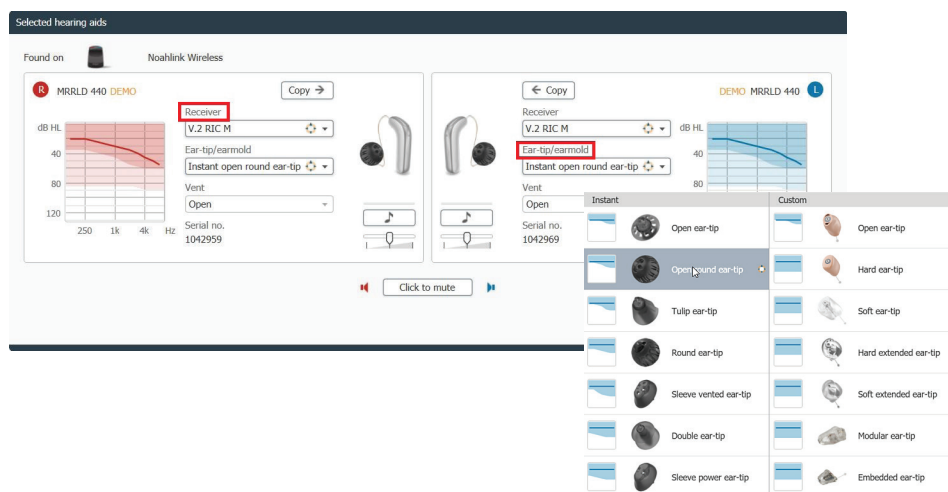

COMPASS GPS FITTING QUICK GUIDE

WIDEX HEARING AIDS

4. When aids are selected and connected for the first time, you will be prompted to select the experience level of the hearing aid wearer.

This allows you to meet the hearing aid wearer where they are with their hearing needs.

TIP: If probe microphone measures will be done, select **Experienced**.

Hearing aid experience

Select the level of experience your client has with hearing aids

The gain settings will be optimized according to your client's experience level.

Select **Inexperienced** if your client is a first-time user or a short-term user.
Select **Experienced** if your client is a long-term user.

Inexperienced


Experienced

Cancel

Feedback Test

5. The Widex Feedback Test ensures the wearer's comfort while still optimizing speech audibility and preserving soft sounds. Wearers need to sit quietly during the test. Press the **start button** to begin.

Note: The feedback test should be run for all Widex fittings.

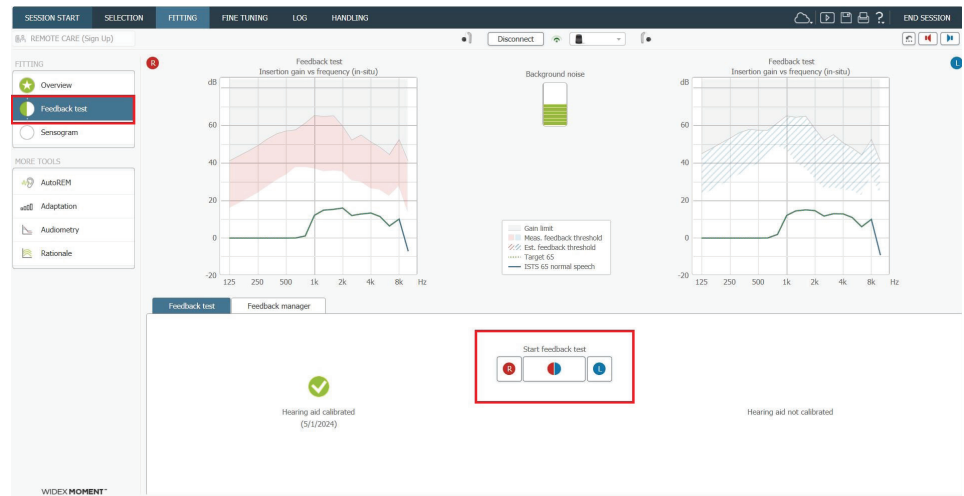


Sensogram

Recommended for the most optimized fitting and to minimize fine tuning.

6. The Sensogram provides a precise in-situ testing option. To perform a Sensogram, simply **highlight each frequency** and present the tone as you would for a hearing test. Monitor room noise to maintain testing accuracy.

TIP: For best results, do not override the measured response. Leave the colored bar at the patient's measured in-situ threshold.

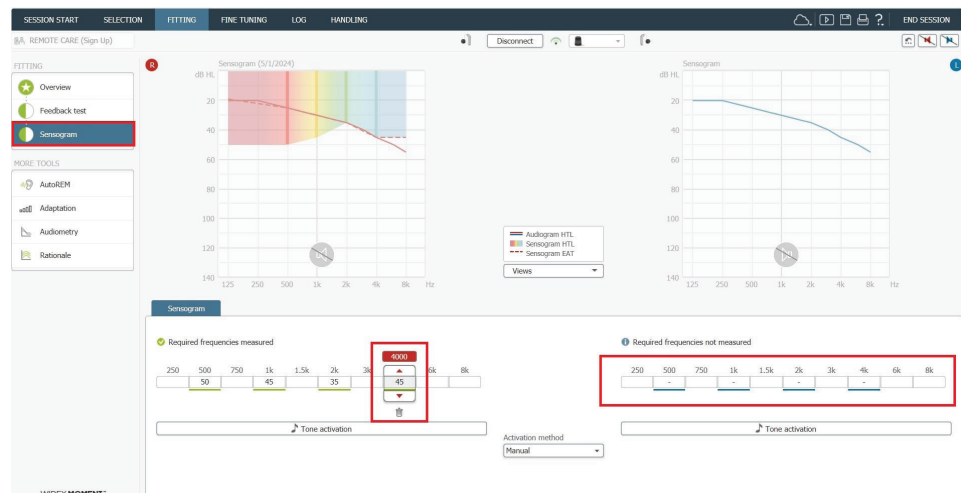


COMPASS GPS FITTING QUICK GUIDE

WIDEX HEARING AIDS

Adaptation

7. The Widex Adaptation manager has four available levels. Level 4 is equal to the full prescribed gain (adaptation off). To address patient comfort, select a level from 1-3 and choose if the hearing should automatically increase over time.


TIP: Most new hearing aid wearers prefer level 2.

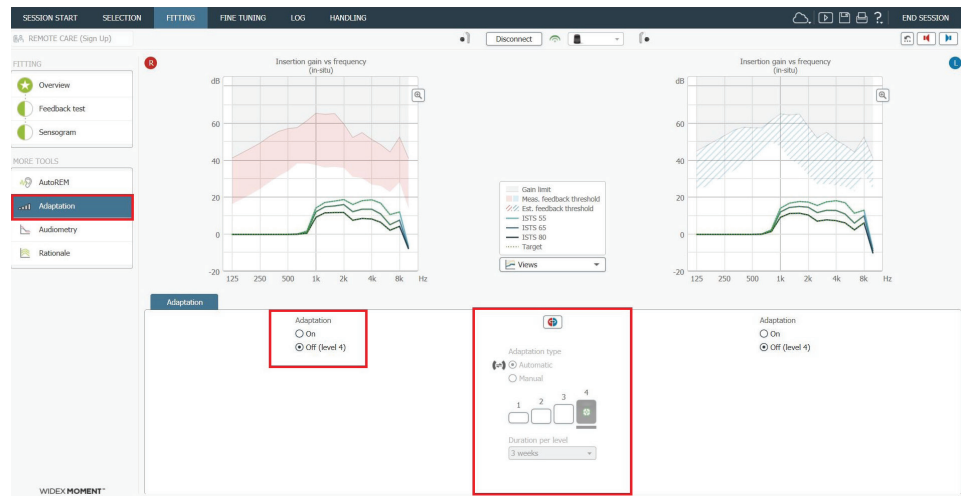

Program Manager


8. Widex offers two possible automatic program options. All Widex fittings have a Universal program. A PureSound program will be automatically added if the wearer's audiogram is within the recommended fitting range. Click on each program icon to demonstrate the sound options to the wearer. Put the wearer's preferred automatic program in the P1 slot in Program manager.

TIP: A PureSound program can be added manually. This option can be considered for newer users that may prefer PureSound but fall outside the recommended (mild/moderate) range.

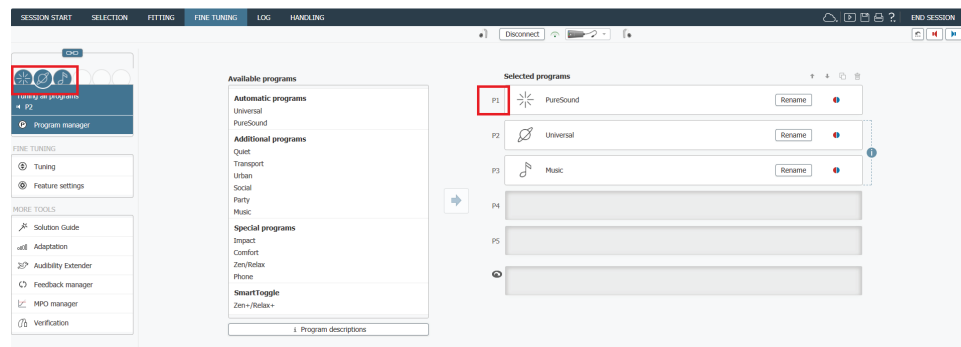

End Session

9. Review the Handling tab for any desired changes to acoustic/light indicators and user controls. Options will vary by model of hearing aid.
10. Click **End Session** to complete your fitting and save the fitting to Noah.

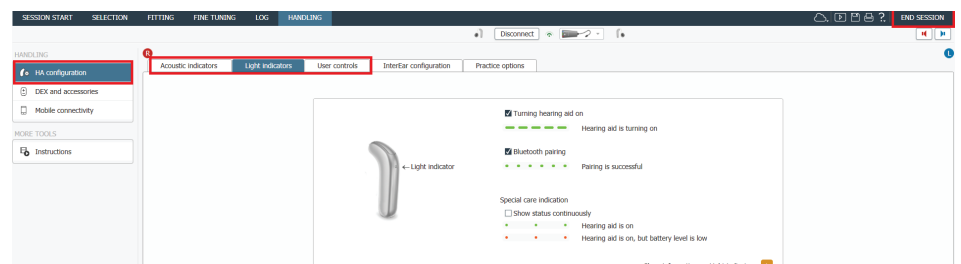


COMPASS GPS FITTING QUICK GUIDE

WIDEX HEARING AIDS

Helpful Fine Tuning Tips


Compass GPS offers three flexible ways to fine tune:

- Sound classes** that can adjust automatic programs for audibility and comfort without adjusting overall gain.
- Customized programs** that can be added for additional support in certain situations.
- Supergain Fast** can be used for increased feedback control.

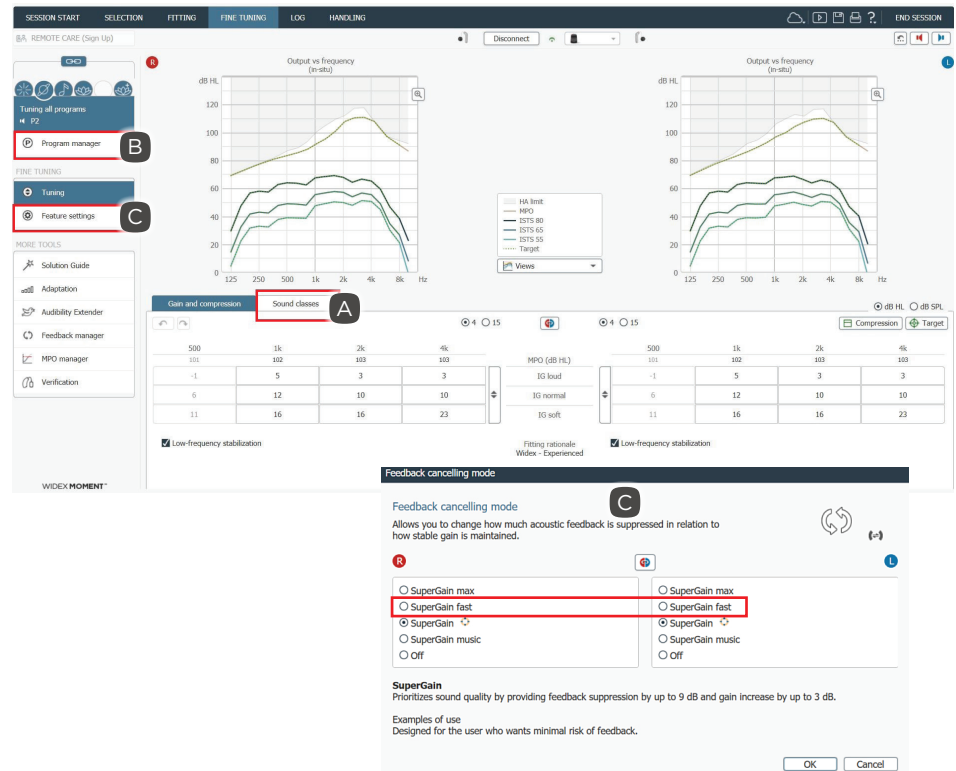

Use the Log tab to personalize both counseling and follow up care.

- Sound Diary** can note manual changes made by the wearer and provides the ability to adjust the sound classes accordingly.
- Real Life Insights** allows viewing of favorite programs in the app if the wearer has created any.

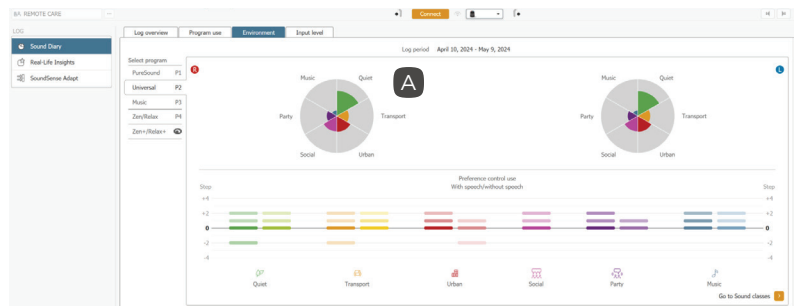

Click **End Session** after finalizing programming changes.

