



1

Learning Outcomes

- 1 Identify 3 feature differences between Beltone Boost Max S and our legacy, Beltone Boost Ultra product.
- 2 List the wireless connectivity options for Beltone Boost Max S.
- 3 List 3 updates to our Beltone Solus Max 2.2 fitting software.



2


Beltone Boost Max S (BMS 95-DWHC)



Tech levels 17, 9, 6 and 4

3


Our smallest super power

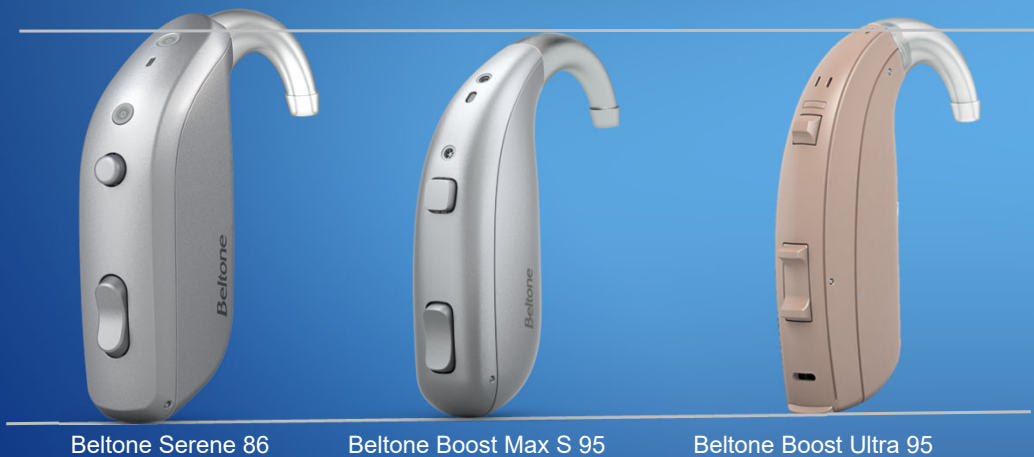

4

Fitting range


Metal Hook - Standard

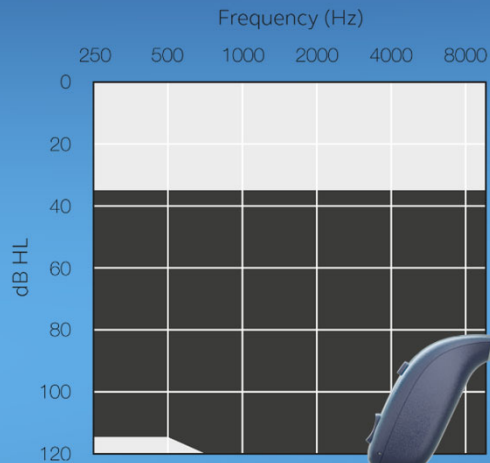

Belton
UNIVERSITY

84 - 140 dB gain

Plastic Hook - Optional

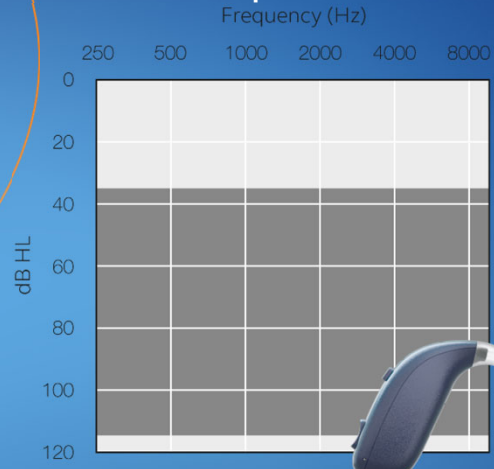

79 - 140 dB gain

5

Exchangeable microphone filters



Belton
UNIVERSITY

6

Premium colors



All-weatherproof IP68-rated and protective nanocoating in and out




7

Boost power



Hours of use

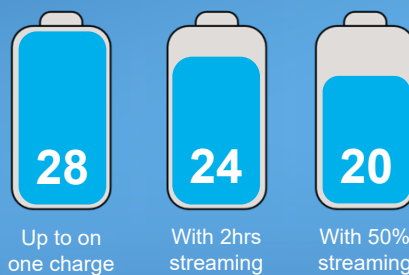

Fast Charging

4 hour = full charge
1 hour = ~7 hours of use
10 min = ~1 hours of use




GN Proprietary data on file (2025)

8

Boost hearing in noise

To help your patients thrive
in the most challenging
listening environments

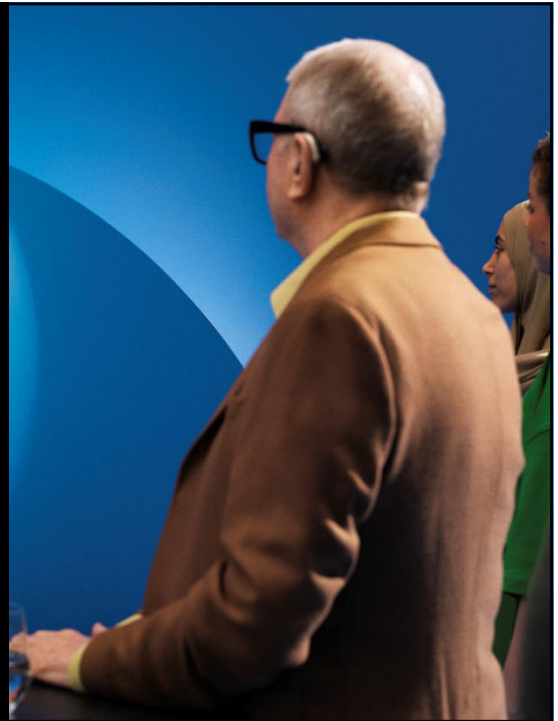

9


CrossLink Directionality 4

Industry best at providing
access to surrounding sound¹

90%



Preference for listening in
CrossLink Directionality 4
over legacy²

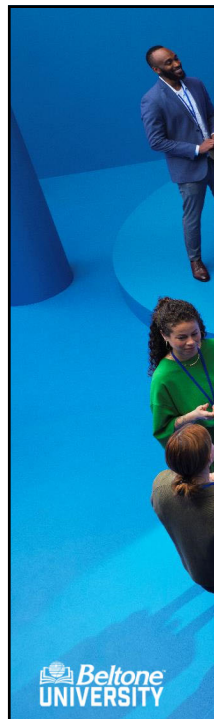

 **Beltone**
UNIVERSITY

¹ Groth et al (2023) + Jespersen et al (2024). ² Jespersen & Skovlund (2025)

10

Clear Focus

Industry's most narrow beamformer¹

2.7dB


SNR improvement for people with severe to profound hearing loss compared to legacy²

8.3dB

SNR improvement for people with severe to profound hearing loss compared to omnidirectional²

93%

Preference for hearing in noise for those with severe to profound hearing loss²

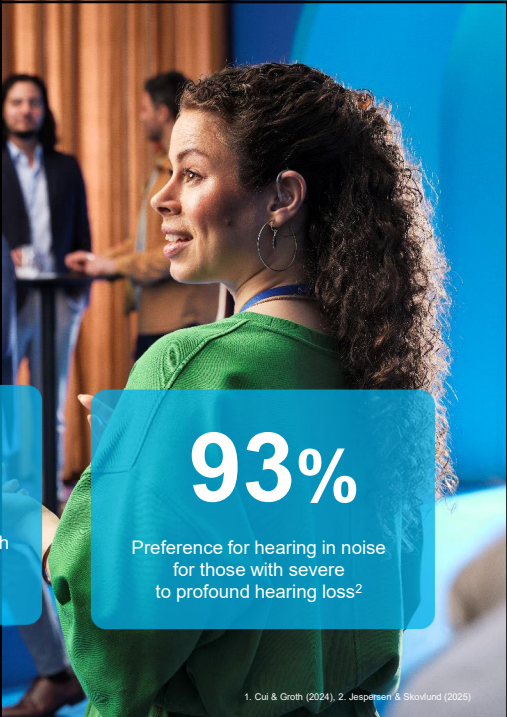

Belton
UNIVERSITY

1. Cui & Groth (2024), 2. Jespersen & Skovlund (2025)

11

Environmental Classifier



Uses AI to scan the environment in real time and classify sound into seven distinct categories



Activates appropriate features based on the amount of noise and the presence of speech



Helps users to hear their best automatically

91%

Accuracy in Environmental classification¹



Belton
UNIVERSITY

1. Groth et al. (2023)

12

Feedback Eraser 2

World's best feedback management¹

Added stable gain (dB)

BELTONE	BRAND A	BRAND B	BRAND C	BRAND D
26	20	17	17	15

Belton
UNIVERSITY

Up to **6dB**

More added stable gain compared to other leading super power solutions¹


1. Cui & Groth (2025)

13

Beltone Boost Max S Technology Differences					
Tech level		17	9	6	4
Hear all around	Spatial awareness	CrossLink Directionality 4	CrossLink Directionality 3	CrossLink Directionality 3	CrossLink Directionality 2
	Hear in noise	Clear Focus	Ultra Focus 2	Ultra Focus	Smart Beam
Personalization	1:1 conversations	Smart Beam Steering	Smart Beam Steering	Smart Beam Steering	Smart Beam Steering
	Group conversations				
Comfort	Number of channels	17	14	12	12
	Adjustable Mixing Point Frequency	✓	✓	✓	✓
	Environmental noise reduction	Per Environment	Strong	Considerable	Moderate
	Impulse noise reduction	3 settings	3 settings	2 settings	2 settings
	Feedback Cancellation	Whistle Stop	Feedback Eraser 2	Feedback Eraser 2	Feedback Eraser 2

14

Boost connectivity

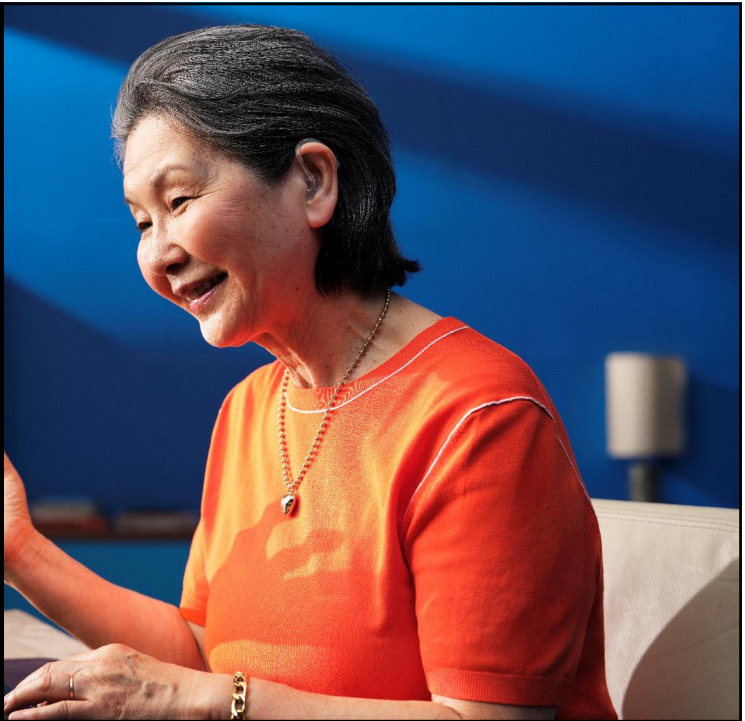

15

Boost entertainment with the TV-Streamer+



Easy pairing with OneTouch Pairing




AutoStreaming makes streaming effortless



Adjust the TV volume and surrounding sounds with the Beltone HearMax app



Supports Auracast

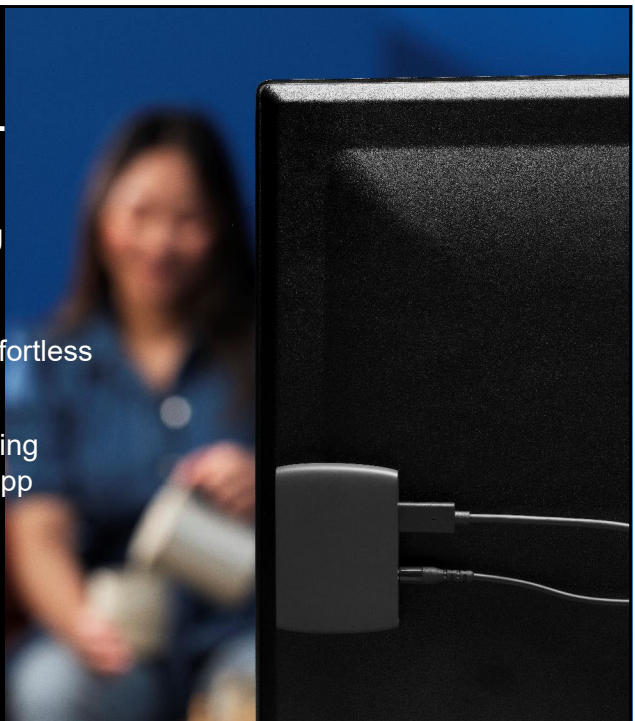

16

Boost conversations with a Multi-Mic+

Capture every word up to 25m away*

9.9dB

SNR improvement compared to omni for Super Power hearing aid users¹

+6dB

Additional SNR improvement than with Clear Focus for Super Power hearing aid users¹



1. For people with severe to profound hearing loss. Jespersen & Rubachandran (2025)


17

Boost streaming using the Auracast Assistant

World's first Auracast Assistant for hearing aids¹


How it works in the Beltone HearMax app:

- 1 Click on "Find a broadcast" at the bottom of the screen
- 2 Choose an Auracast broadcast from the list
- 3 Enjoy sound directly streamed to your hearing aids

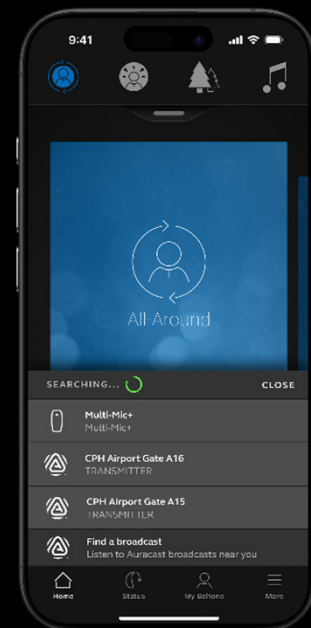

1. GN Proprietary data on file (2025)

18

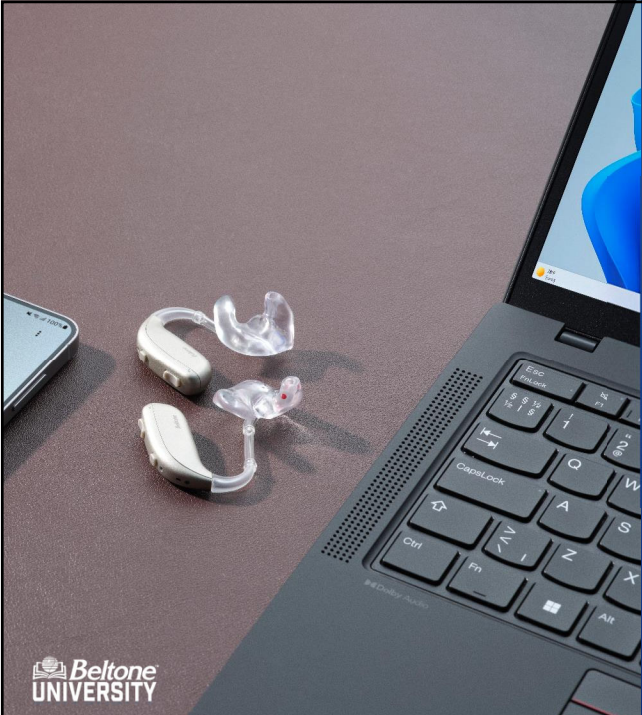

Boost convenience with MultiConnect




Pair Beltone Boost Max S with two smart devices simultaneously for uninterrupted streaming – such as taking a phone call while streaming music or movie audio from a laptop.

19

What's new in Beltone Solus Max 2.2?

1 New discovery screen

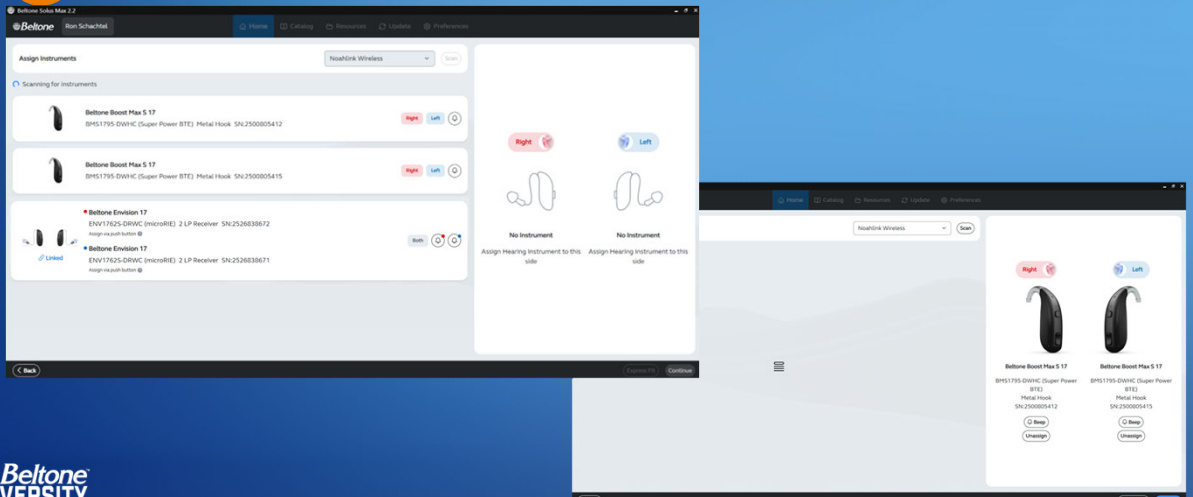


20

2

Updated features listing

me Catalog Resources Update Preferences

Features Specification

Speech

CrossLink Directionality 4

True Focus

Personal Sound ID

Mixing Point Frequency

Beamwidth

Sound Shifter

Comfort

True Noise Reduction

Sound Cleaner Pro

Silencer

Feedback Eraser 2

Wind Noise Reduction

Impulse Noise Reduction

Low Frequency Boost

Automatic Gain

Smart Gain Pro

SATISFY

Linear Compression

Features Specification

Speech

CrossLink Directionality 4

True Focus

Clear Focus

Personal Sound ID

Mixing Point Frequency

Beamwidth

Sound Shifter

Comfort

True Noise Reduction

Sound Cleaner Pro

Silencer

Feedback Eraser 2

Wind Noise Reduction

Impulse Noise Reduction

Low Frequency Boost

Automatic Gain

Smart Gain Pro

SATISFY

21

3

Faster firmware updates

Beltono Solus Max 2.2

Firmware Updates

Right Instrument

Beltono Envision 17 ENV1762S-DRWC (microRIG) 2 LP Receiver

Left Instrument

Beltono Envision 17 ENV1762S-DRWC (microRIG) 2 LP Receiver

Step 2 of 7

Updating firmware for both sides.

The update may take several minutes.


Do not interrupt the update once it has begun.

Save fitting to apply update to the hearing aid(s) before exiting.

22


11

4 New feature for CROS

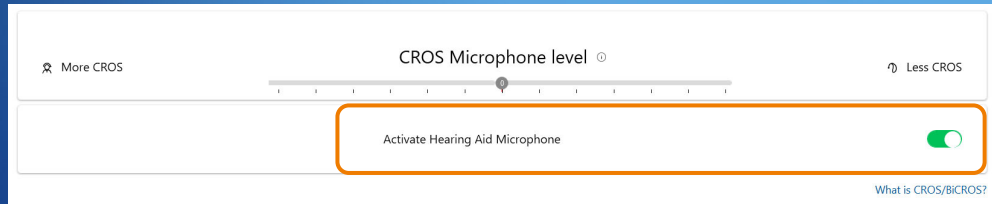

23

5 New accessory firmware

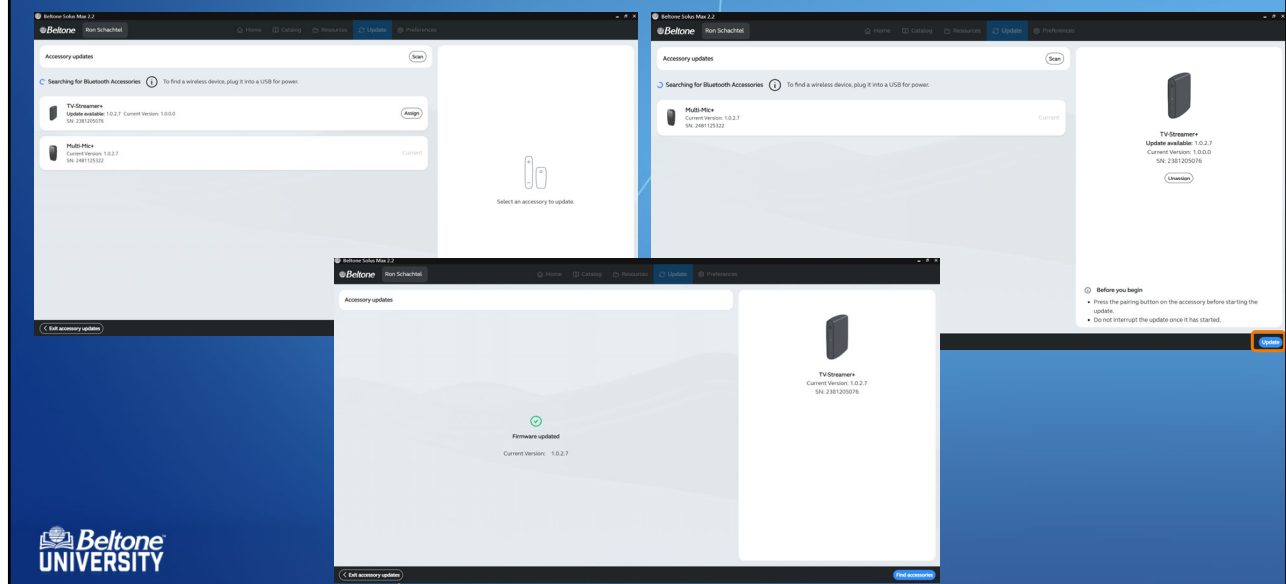

24

6 New test device



Hearing aids remain connected to software throughout the test



Ability to mute hearing aids from the software




Now with directionality tests



Viewable Specs




Hearing aids go back to patient's settings

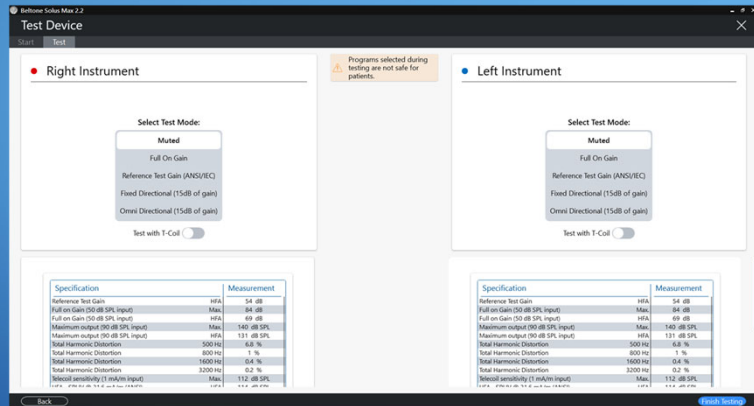

25

Beltone HearMax app: True Noise Reduction control

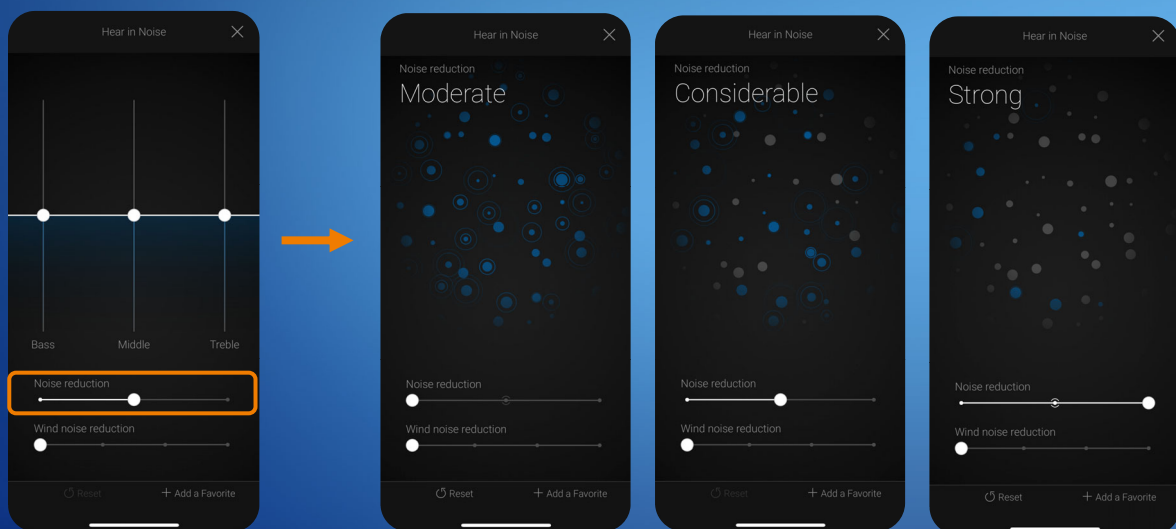

26

Beltone Boost Max S made for hearing in noise



World's smallest rechargeable super power hearing aid¹



Fast charging



CrossLink Directionality 4



28 hrs
All-day battery



World's first super power hearing aids with Bluetooth LE Audio & Auracast⁴



All-weatherproof
IP68 rated



AURACAST

93%

preference for hearing in noise over legacy²



Clear Focus

New premium colors



Exchangeable mic filters



World's best feedback manager³





Multi-Mic+
TV-Streamer+



1. Compared to competing Super Power solutions. GN proprietary data on file (2025). 2. Groth & Cui (2025). 3. Jespersen & Skovlund (2025). 4. GN proprietary data on file (2025)

27

Summary

1

Beltone Boost Max S now runs on our 360 chip. Our 17 level offer the patients access to CrossLink Directionality 4, Clear Focus, and up to 6 dB of stable gain with Feedback Eraser 2 than others.

2

Beltone Boost Max S use Bluetooth LE Audio and can connect to our streaming accessories and Auracast broadcasts.

3

Several screens have new looks in Beltone Solus Max 2.2; the discovery screen, accessory update and catalog list. Also hearing aid firmware updates are fast as they can be completed simultaneously.



28