

Custom Fit, Max Impact: Beltone Commence Customs & Solus Max Updates



Learning Outcomes

1

Explain the hardware changes with the Beltone Commence custom styles.

2


List 3 reasons Beltone Commence is superior to Beltone Rely.

3

Recall the enhancements in Beltone Solus Max 2.3.

CIC Options

- Receiver options

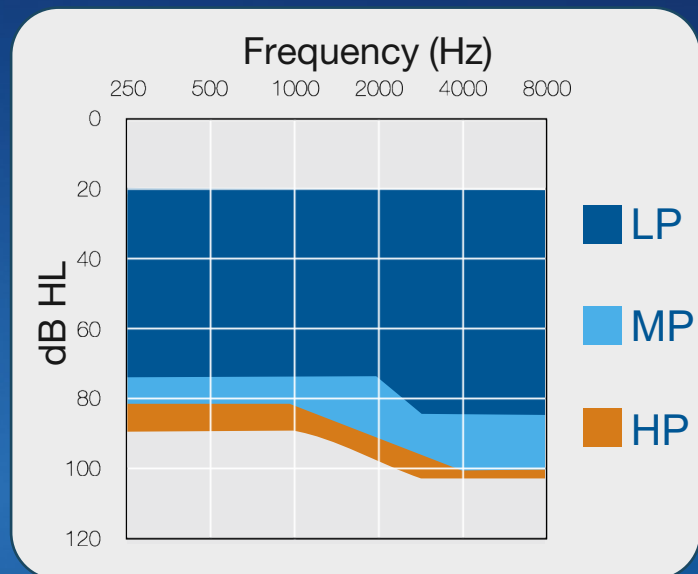

- Optional push button
- 10A zinc air battery
- Cabled programming

- Colors:



Anthracite

Light
Beige

Beige

Medium
Brown


Brown

Dark
Brown


Blue/Red/Clear
shells

ITC Options

- Receivers

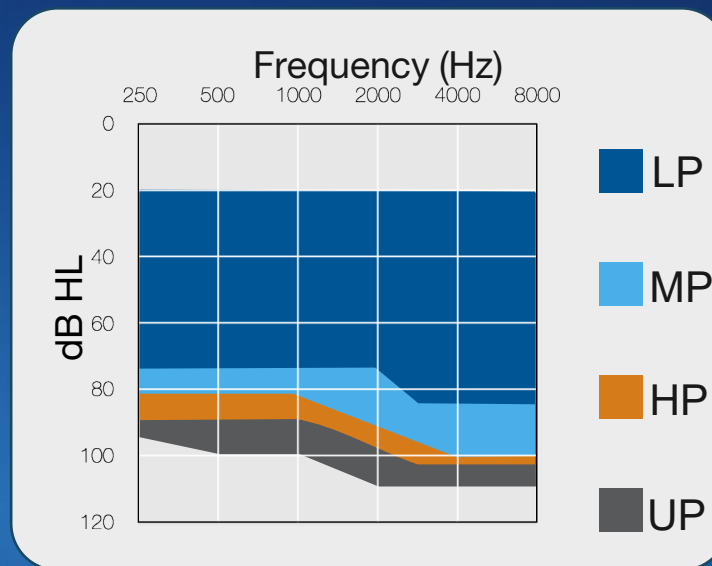

- Push button and optional VC
- 312 battery in shell options:

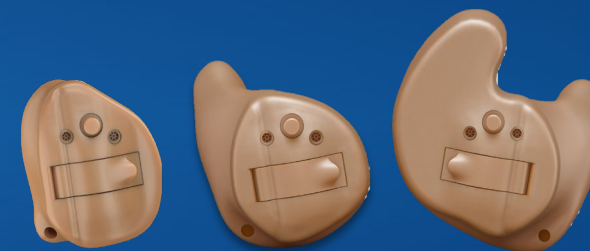

ITC Half shell Full shell

- Colors:



Anthracite Light Beige Beige Medium Brown Brown Dark Brown Blue/Red/Clear shells

MEMS Microphone Matters



More consistent
performance



Better resistance
to moisture and
debris



Compact and
precise

Microphone Filter

Old



New



Beltone Commence RIEs



Wireless Accessories:

Giving patients more control

- Easy OneTouch Pairing
- AutoStreaming
- Beltone HearMax app adjustability
- Supports Auracast™



Multi-Mic+



TV-Streamer+

Bluetooth® LE Audio and Auracast All-inclusive connectivity

**+2.8
million**

ears helped with Auracast
enabled hearing aids
since 2023²

76%

Prefer hearing with
Auracast compared
to teleloop³

+5

New Auracast enabled
venues pop up every
day globally¹



Beltone Commence Connectivity Options

					
	CIC	ITC	62S miniRIE	63 RIE	64 RIE
MFi/Android streaming		●	●	●	●
Hands-free calls for iPhone & iPad*			●	●	●
Bluetooth LE Audio / Auracast		●	●	●	●
Wireless accessories (LE Audio/Auracast)		●	●	●	●
Wireless accessories (2.4 GHz)		●	●	●	●
Connection to software and smartphone		●	●	●	●
Telecoil					●
Automatic receiver detection			●	●	●
Sound Shifter	●	●	●	●	●

Features: Beltone Commence vs. Beltone Rely

Beltone Rely 4

You get the Beltone Commence 3 minus:

- Impulse Noise Reduction
- Smart Beam directionality
- Adjustable Mixing Point Frequency

Beltone Rely 3

You get **less than** Beltone Commence 2 experience,

- Less channels
- No E2E synchronisation

Beltone Rely 2

You get **less than** Beltone Commence 2 experience:

- Less channels
- No streaming to Smartphones



Exchangeable
mic filters*



Bluetooth I.E. Audio,
Auracast &
hands-free calls**



simultaneous
connection
Solus Max &
Smartphone

Beltone Commence 3

- + 12 channels
- + Impulse Noise Reduction
- + Smart Beam directionality
- + Adjustable Mixing Point Frequency
- + Streaming Bass Boost
- + Low Frequency Boost



Beltone Commence 2

- + 10 channels
- + Streaming Bass Boost
- + Low Frequency Boost



Connectivity: Beltone Commence vs. Beltone Rely

Beltone Commence



Beltone Rely



	CIC	ITC	miniRIE	63 RIE	64 RIE	CIC	ITC	ITE	RIE-C	63 RIE	64 RIE	
MFfi/Android streaming		●	●	●	●		●	●	●	●	●	●
Hands-free calls for iPhone & iPad*			●	●	●							
Bluetooth LE Audio / Auracast		●	●	●	●							
Wireless accessories (LE Audio/Auracast)		●	●	●	●							
Wireless accessories (2.4 GHz)		●	●	●	●	●	●	●	●	●	●	●
Connection to software & smartphone		●	●	●	●							
Telecoil					●					●	●	●
Automatic receiver detection (RIE)			●	●	●							
Sound Shifter	●	●	●	●	●							

Firmware update

Right Instrument

Beltone Envision 17 ENV1762S-DRWC (microRIE) 2 LP Receiver



microRIE

Serial No: 2526838672

Firmware version: 10.60.1.0

Right Instrument

The update may take several minutes.
Do not interrupt the update once it has begun.
Ensure your rechargeable hearing aids have sufficient battery charge, and are removed from the charger.
Insert fresh battery in non-rechargeable hearing aids.
Hearing instruments must not be worn on the ears during firmware update process.

Left Instrument

Beltone Envision 17 ENV1762S-DRWC (microRIE) 2 LP Receiver



microRIE

Serial No: 2526838671

Firmware version: 10.21.1.1

A newer version of firmware is available for this instrument:
Firmware version: 10.60.1.0

Firmware Update Available

☐ Update Firmware

☒ Continue Without Updating Firmware

The update may take several minutes.
Do not interrupt the update once it has begun.
Ensure your rechargeable hearing aids have sufficient battery charge, and are removed from the charger.
Insert fresh battery in non-rechargeable hearing aids.
Hearing instruments must not be worn on the ears during firmware update process.

Back

Continue

Right Instrument

on 17 ENV1762S-DRWC (microRIE) 2 LP Receiver



microRIE

Serial No: 2526838672

Firmware version: 10.21.1.1

A newer version of firmware is available for this instrument:
Firmware version: 10.60.1.0

Firmware Update Available

☒ Update Firmware

☐ Continue Without Updating Firmware

The update may take several minutes.
Do not interrupt the update once it has begun.
Ensure your rechargeable hearing aids have sufficient battery charge, and are removed from the charger.
Insert fresh battery in non-rechargeable hearing aids.
Hearing instruments must not be worn on the ears during firmware update process.

Left Instrument

Beltone Envision 17 ENV1762S-DRWC (microRIE) 2 LP Receiver



microRIE

Serial No: 2526838671

Firmware version: 10.21.1.1

A newer version of firmware is available for this instrument:
Firmware version: 10.60.1.0

Back

Continue

Firmware update

GN Online Services Login    

 Menu 

Patient

Speech Scores

RECD | REUG

Instruments

Express Fit

Reset to Initial Fit

Reconfigure 

Test Device

Fitting

Beltone Solus Max Guide

Notifications

Firmware Updates

Other

Sound Player


Audiogram History

Beltone Solus Max 2.3

Audiogram Mismatch

● Right Instrument

Beltone Envision 17 ENV1762S-DRWC (microRIE) **2 LP** Receiver

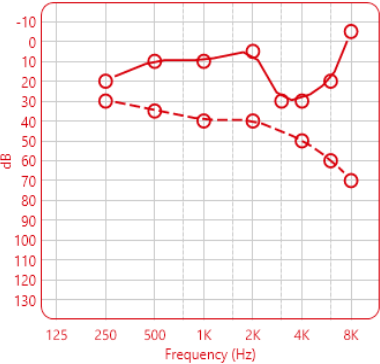

⚠ Audiogram Mismatch

—● Audiogram History

5/22/2023 5:04 PM

○ Use Instrument Audiogram

5/22/2023 5:04 PM

9/15/2025 6:09:45 PM


12/13/2016 1:37 PM

5/31/2016 1:42 PM

1/11/2016 10:23 AM

● Left Instrument

Beltone Envision 17 ENV1762S-DRWC (microRIE) **2 LP** Receiver

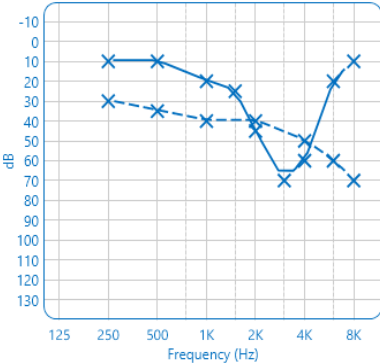

⚠ Audiogram Mismatch

—● Audiogram History

5/22/2023 5:04 PM

○ Use Instrument Audiogram

5/22/2023 5:04 PM

9/15/2025 6:09:45 PM

Back

Continue

Audiogram History

Patient
Instrument Set-Up
Gain Adjustments
Advanced Features
Device Controls
Finalize

Profile
In-Situ Fitting
Data Logging
Remote Care

Prescribed First Fit ⓘ

Target Rule
BAFA

? Audiogram 5/22/2023

Change Audiogram

Change Audiogram

Choose the audiograms to use for this fitting.

● Right Instrument

Beltone Envision 17 ENV1762S-DRWC (microRIE) 2 LP Receiver

Audiogram History 5/31/2016 1:42:03 PM

Use Instrument Audiogram 5/22/2023 5:04:18 PM

9/15/2025 6:09:45 PM

12/13/2016 1:37:39 PM

5/31/2016 1:42:03 PM

1/11/2016 10:23:38 AM

● Left Instrument

Beltone Envision 17 ENV1762S-DRWC (microRIE) 2 LP Receiver

Audiogram History 5/22/2023 5:04:18 PM

Use Instrument Audiogram 9/15/2025 6:09:45 PM

Cancel
Continue

Programming Interface

Beltone Solus Max 2.3

Bob Dome

Home Catalog Resources Update Preferences

Assign Instruments

Noahlink Wireless

Stop

Scanning for instruments

If the patient's hearing aids cannot be found, please notify the patient to turn off the Bluetooth on mobile phones that are paired to the hearing aids.

Noahlink Wireless

Noahlink Wireless

Speedlink

NOAHlink

HI-PRO 2

Scan

Speedlink



Scan

Back

Express Fit

Continue

Software Version Number

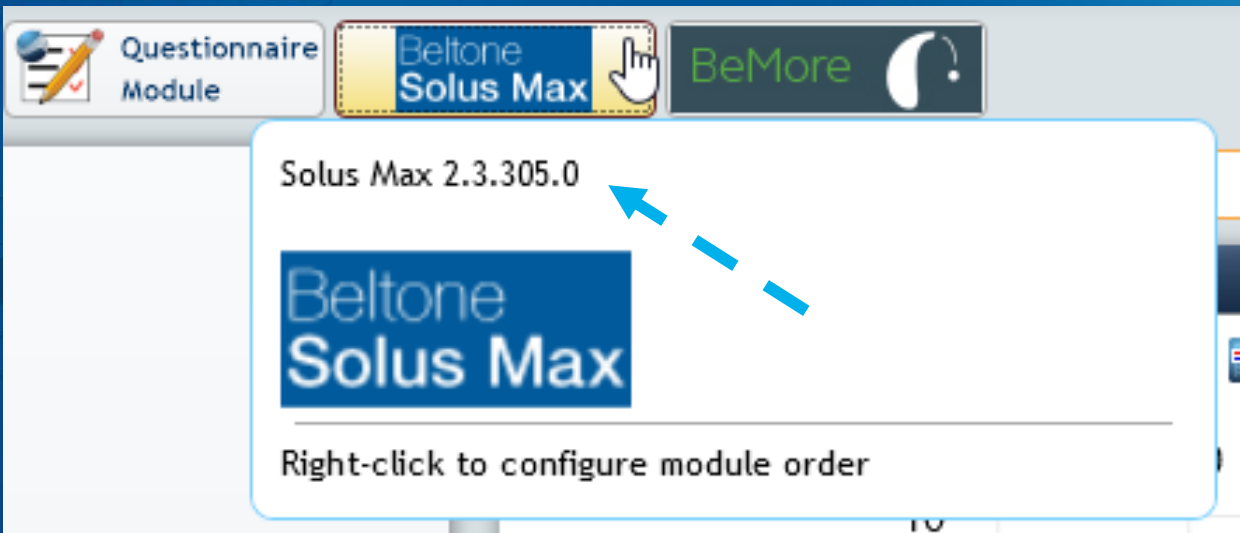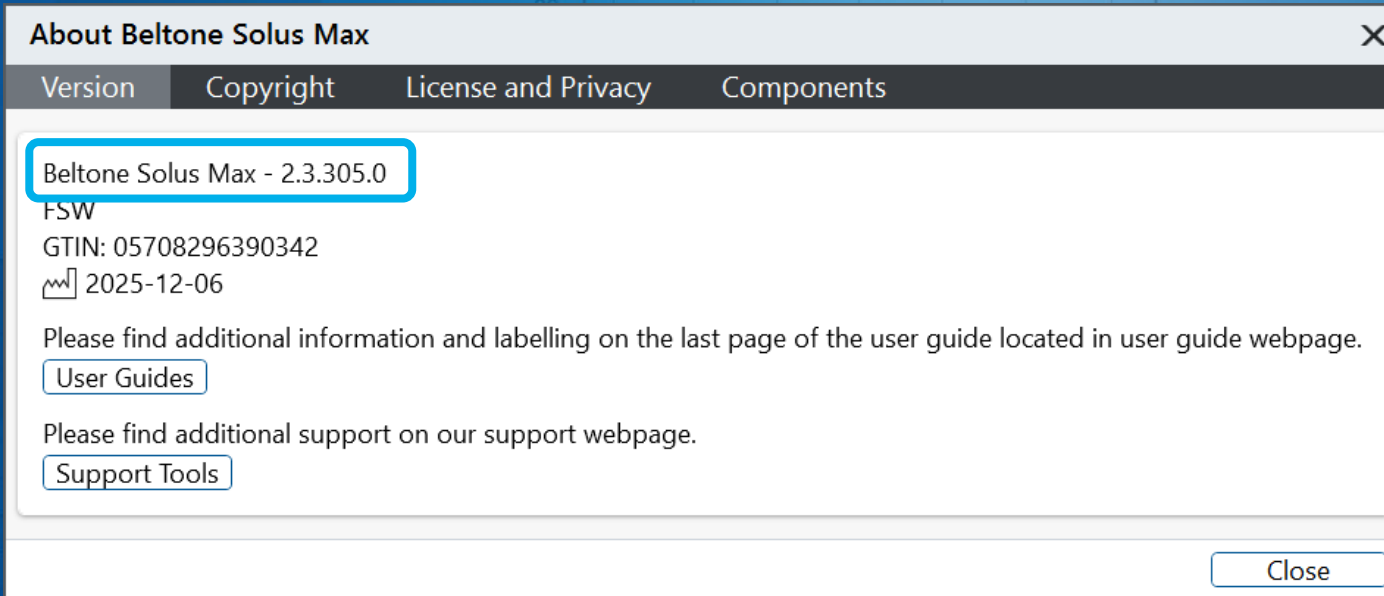


Beltone Commence



**Hearing clarity and
connectivity
for all**

Offered in RIE &
custom styles

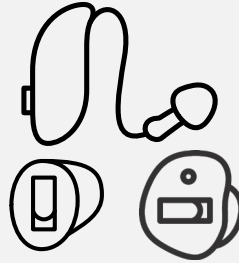

Best feedback
manager⁵



24 hrs

All-day battery



World's
first
Auracast
Assistant⁴



All-weatherproof
IP68 rated



Compatible with new AllFit domes



BEST-IN-CLASS

Hearing quality¹
Hearing comfort²
Bluetooth connectivity³



Impulse Noise
Reduction

Color options



Replaceable
microphone
filters, RIE &
ITC



Enjoy the
Multi-Mic+
TV-Streamer+

Learning Outcomes

1

Beltone Commence custom styles new use MEMS microphones.

2

Beltone Commence uses our 360 chip, has more advanced features, uses BTLEA and Auracast and Hands-free calls

3

Beltone Solus Max 2.3 allows for audiogram selection, firmware update in the connection flow and they ability to change the programing interface from the connection screen.