

7 reasons why ReSound Smart Fit™ will enrich your fitting experience



Download the latest fitting software on pro.resound.com



As a hearing care professional, you spend a lot of your day in fitting software. You deserve the best experience. ReSound Smart Fit fitting software keeps evolving and improving based on your feedback and user testing to better meet your needs. Enjoy smarter and more intuitive software, every time you fit with the latest ReSound Smart Fit.



Spend less time clicking

Enjoy a simple experience when navigating the home screen and preferences.



Enjoy a more intuitive experience

Simplified main navigation with a vertically organized design for a more intuitive fitting experience.



Improved features

Instrument card is enhanced with a new battery level display, now includes serial number and firmware version, and available on all the screens so that the hearing aid information as well as connection and calibration are readily accessible. Enhancements on Test Device feature were made to include display of battery status and to reduce scrolling.



Complete verification testing

Use directionality presets in our fitting software to verify directional microphone performance in the test box and complete all verification measurements while staying connected to the hearing aids – in 50% less time than before.²

ReSound Smart Fit is the world's first fitting software supporting NAL-NL3™, the next-generation prescriptive fitting system by National Acoustics Laboratory.



Personalized onboarding for your patients

Provide your patients with personalized guidance by printing onboarding reports with all the features and guidance they need.



Revamped home screen

Stay up to date on software, accessory firmware, and all our latest materials, now directly accessible from the home screen.



Adapts to your way of working



Select your favorite settings to become part of your standard fitting flow. For your convenience, ReSound Smart Fit will also keep track of your fittings and learn from them. If it recognizes that you do the same thing consistently, it will automatically suggest that preference be set for your fitting flow moving forward.

¹Compared to ReSound Smart Fit 1.17, technical data on file

²50% time saving compared to ReSound Smart Fit 2.1 (from 5.36 mins to 2.54 mins)

6 easy steps to fit ReSound

- 1 Launch ReSound Smart Fit and put the hearing aids in discovery mode by removing them from the charger or closing the battery door for non-rechargeable models. Click **Connect**.
- 2 The hearing aids will appear in the **Assign Instruments** section. Use the **Push to Assign** feature to assign to the respective right and left ear placement. Click **Continue** in the lower right of the screen. Select appropriate physical properties and click **Continue** to proceed through connection process.
- 3 Calibrate the hearing aids on the patient's ears to activate DFS Ultra III. Click **Continue** in the lower right of the screen. **Note:** Hearing aids will automatically unmute once calibrated.
- 4 Navigate to **Patient > Profile** screen, customize the Prescribed First Fit selections as appropriate:
 - a Select Patient **Experience Level** based on previous experience with amplification and to maximize acceptance of hearing aids.
 - b Navigate to **Gain Adjustments**
- 5 Fine-tune hearing aids using **Fine Tuning** under **Gain Adjustments**:
 - a Add or remove programs from this screen.
 - b Once programming is set, navigate to **Device Controls** to customize the manual push button, demonstrate beeps or pair accessories.
- 6 Click **Save** at bottom right of screen. Devices will reboot and play the start-up melody upon exiting the software.

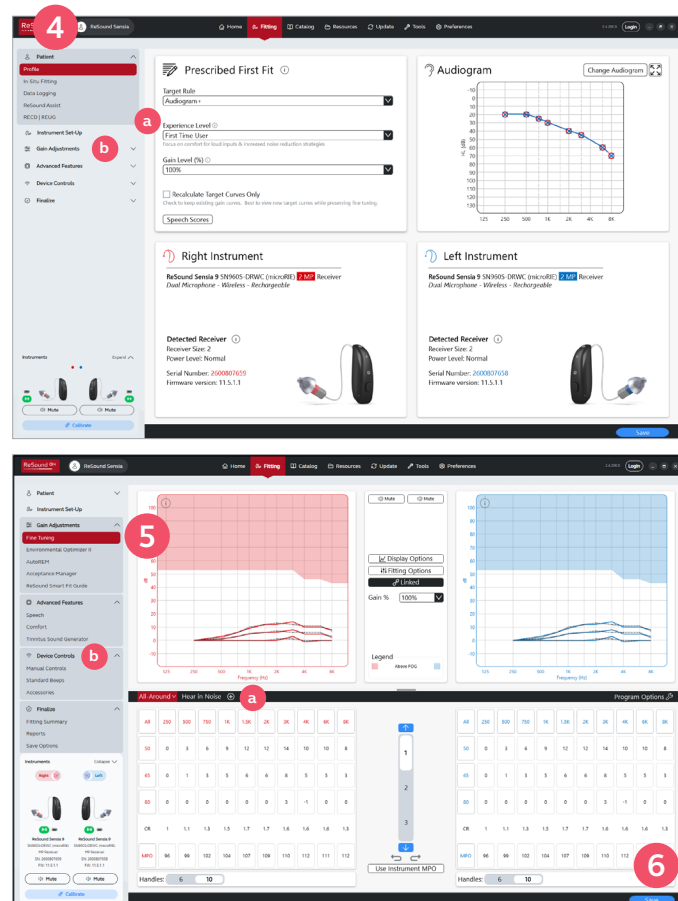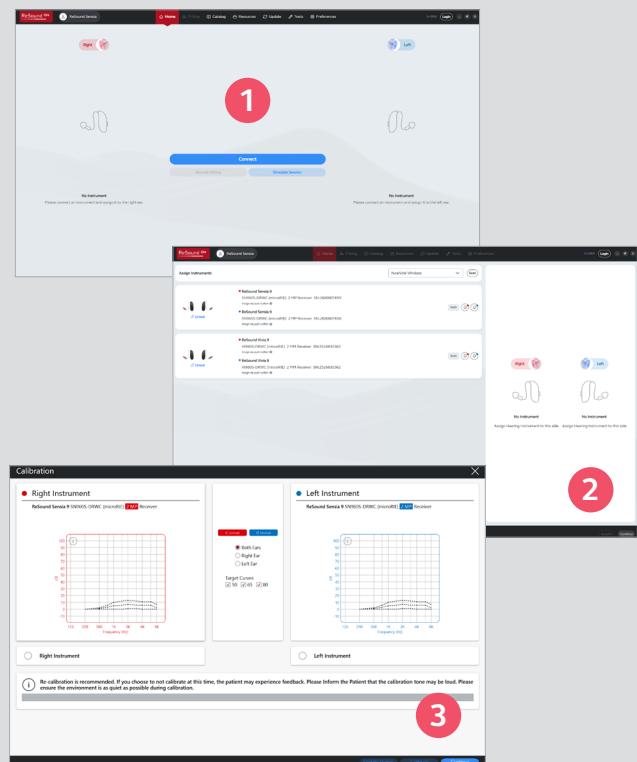

Manufacturer according to FDA:
GN ReSound North America
 5005 Dean Lakes Blvd.
 Shakopee, MN 55379
 USA
 1-800-248-4327
 pro.resound.com

Manufacturer according to Health Canada:
ReSound Canada
 2 East Beaver Creek Road, Building 3
 Richmond Hill, ON L4B 2N3
 Canada
 1-888-737-6863
 pro.resound.com