

What's new: ReSound Smart Fit™ 2.4

ReSound GN




Enjoy a smarter and more intuitive software

We are continuing our journey to make ReSound Smart Fit your preferred software that fits your way of working. Your feedback has been integral to ReSound Smart Fit 2.4 with improvements to the main navigation and enhancements for improved fittings. We continue to make improvements to create an even more intuitive and easy fitting experience for you.

ReSound Smart Fit 2.4 includes:


- Support for new ReSound Sensia™ models
- Support for new ReSound Sensia Bimodal and ReSound Enzo™ IA Bimodal models
- Simplified main navigation with a vertically organized design for a more intuitive fitting experience
- Enhancements on Test Device feature, which includes display of battery status, test standard (IEC or ANSI) and full visibility of the selected standard without scrolling

**Upgrade your fitting software
today by visiting pro.resound.com.**

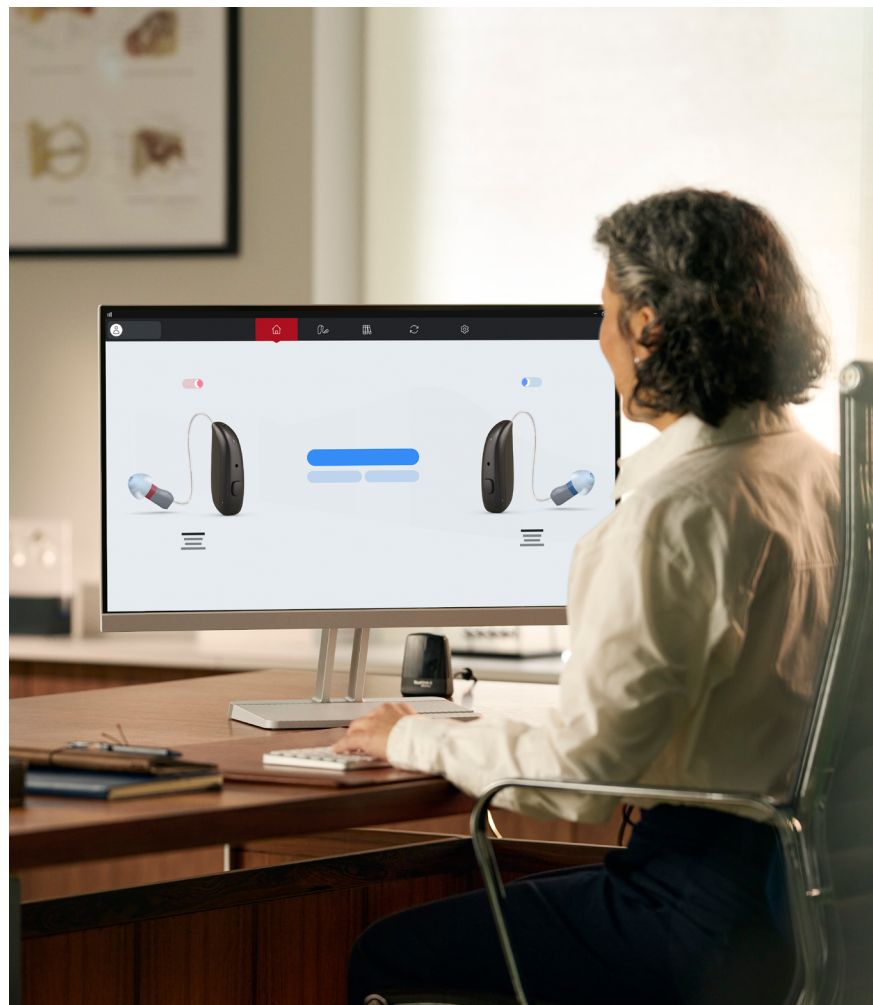


6 easy steps to fit ReSound

- 1 Launch ReSound Smart Fit and put the hearing aids in discovery mode by removing them from the charger or closing the battery door for non-rechargeable models. Click **Connect**.
- 2 The hearing aids will appear in the **Assign Instruments** section. Use the **Push to Assign** feature to assign to the respective right and left ear placement. Click **Continue** in the lower right of the screen. Select appropriate physical properties and click **Continue** to proceed through connection process.
- 3 Calibrate the hearing aids on the patient's ears to activate DFS Ultra III. Click **Continue** in the lower right of the screen. **Note:** Hearing aids will automatically unmute once calibrated.

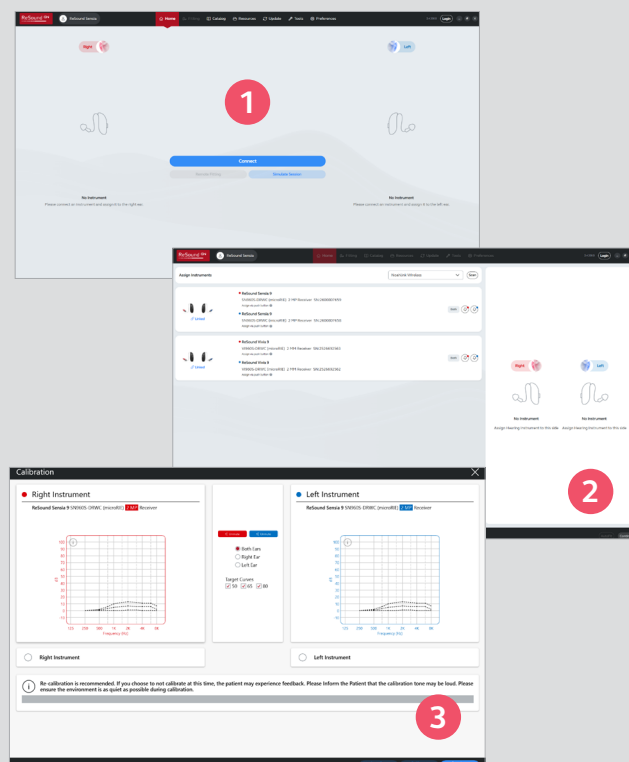

- 4 Navigate to **Patient > Profile** screen, customize the Prescribed First Fit selections as appropriate:
 - a Select Patient **Experience Level** based on previous experience with amplification and to maximize acceptance of hearing aids.
 - b Navigate to **Gain Adjustments**
- 5 Fine-tune hearing aids using **Fine Tuning** under **Gain Adjustments**:
 - a Add or remove programs from this screen.
 - b Once programming is set, navigate to **Device Controls** to customize the manual push button, demonstrate beeps or pair accessories.

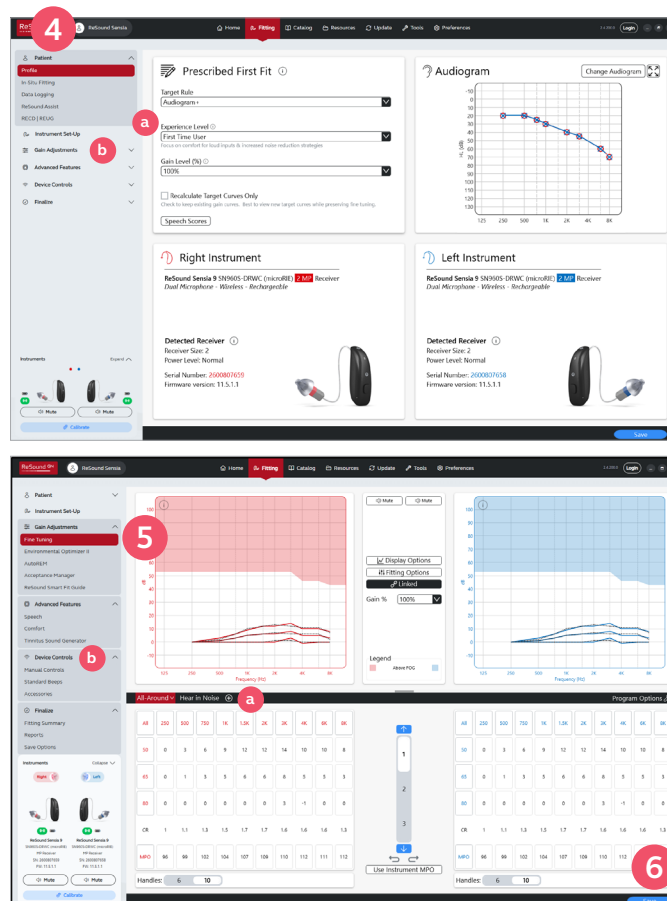

- 6 Click **Save** at bottom right of screen. Devices will reboot and play the start-up melody upon exiting the software.

Manufacturer according to FDA:

GN ReSound North America
5005 Dean Lakes Blvd.
Shakopee, MN 55379
USA
1-800-248-4327
pro.resound.com

Manufacturer according to Health Canada:

ReSound Canada
2 East Beaver Creek Road, Building 3
Richmond Hill, ON L4B 2N3
Canada
1-888-737-6863
pro.resound.com