

Bringing it dome
...because fit happens



Learning Outcomes

1

Describe the improvements seen with our new AllFit dome portfolio.

2

List the different occlusion properties for the different AllFit dome styles.

3

Identify the appropriate AllFit dome based on ear canal anatomy, hearing loss configuration, and occlusion considerations.

AllFit Domes

Improved

**Acoustic properties
& durability**

Thicker silicone for better stability and reduced risk of feedback.

97%

coverage based on thousands of real ear impressions¹

Additional wax protection

with new dome tip

**NEW
Sleeve dome**

Designed for a custom-like fit and optimized for M&RIE



Enhanced low frequency sound quality
(More bass response below 2kHz)




New
Open dome
design

5 sizes




AllFit Open vs SureFit 3 Open

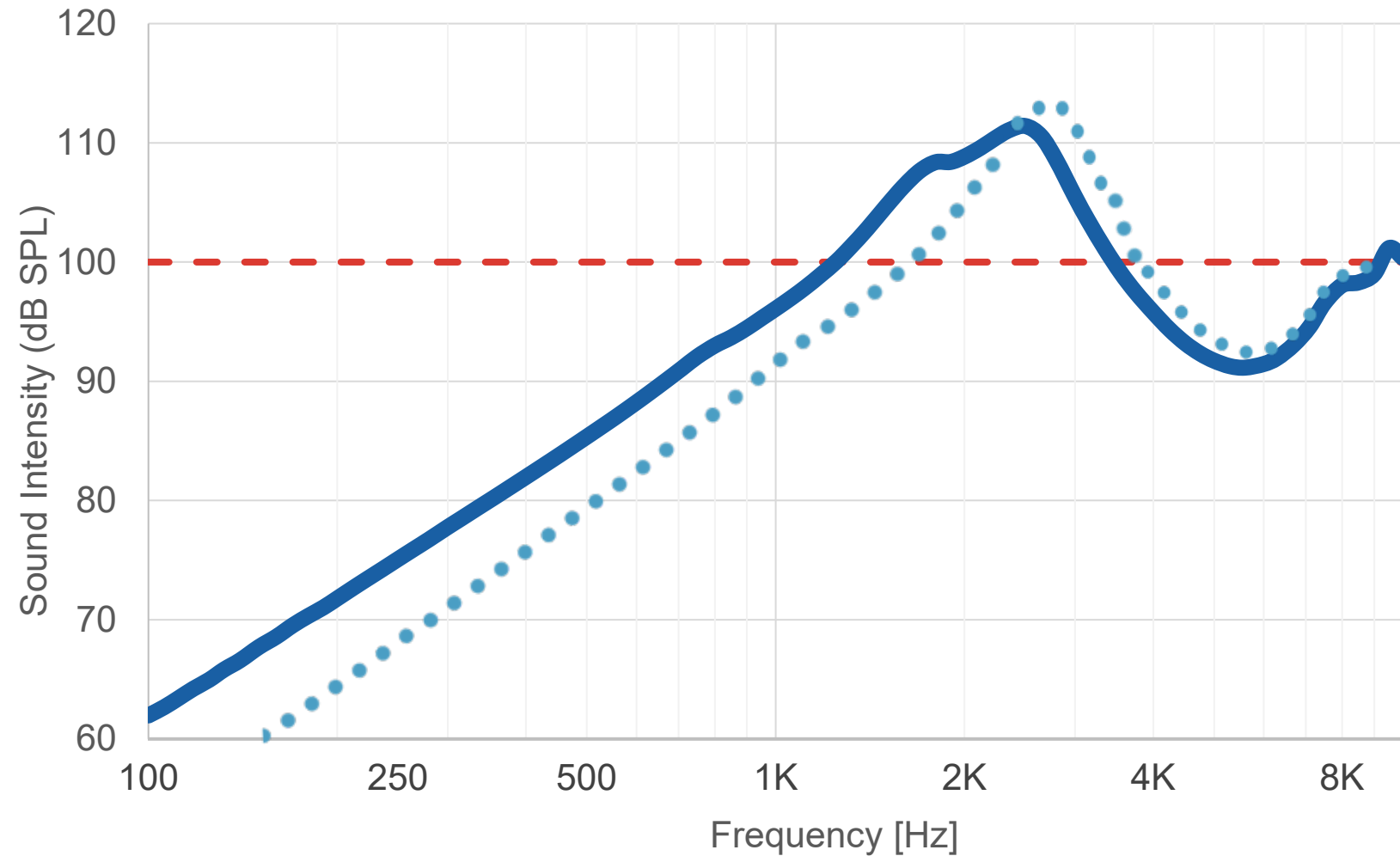


Open Dome

The most vented dome in the collection to provide the greatest relief from occlusion.

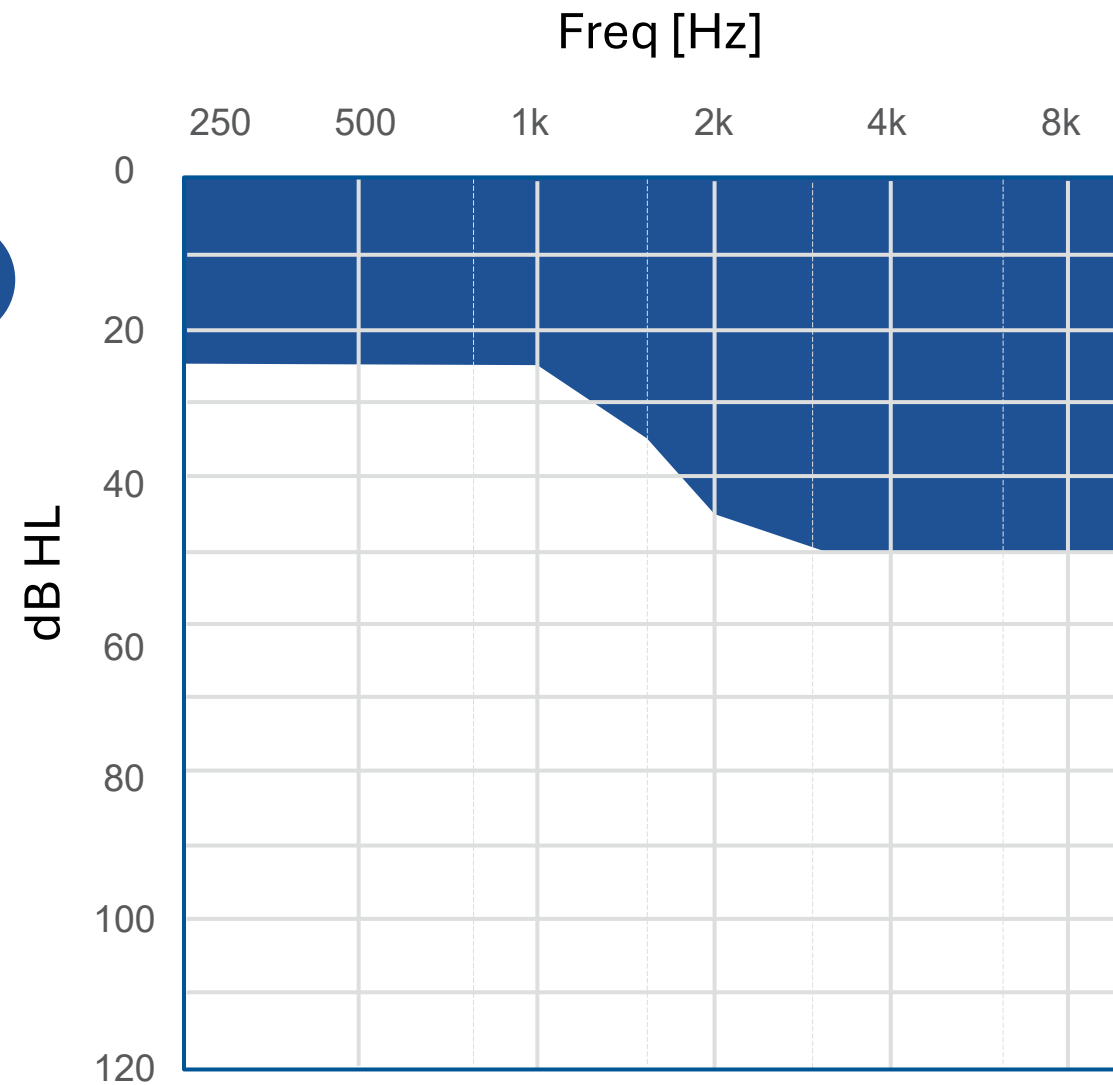


Unique and flexible shape

Good choice for a variety of ear canal shapes, especially narrow

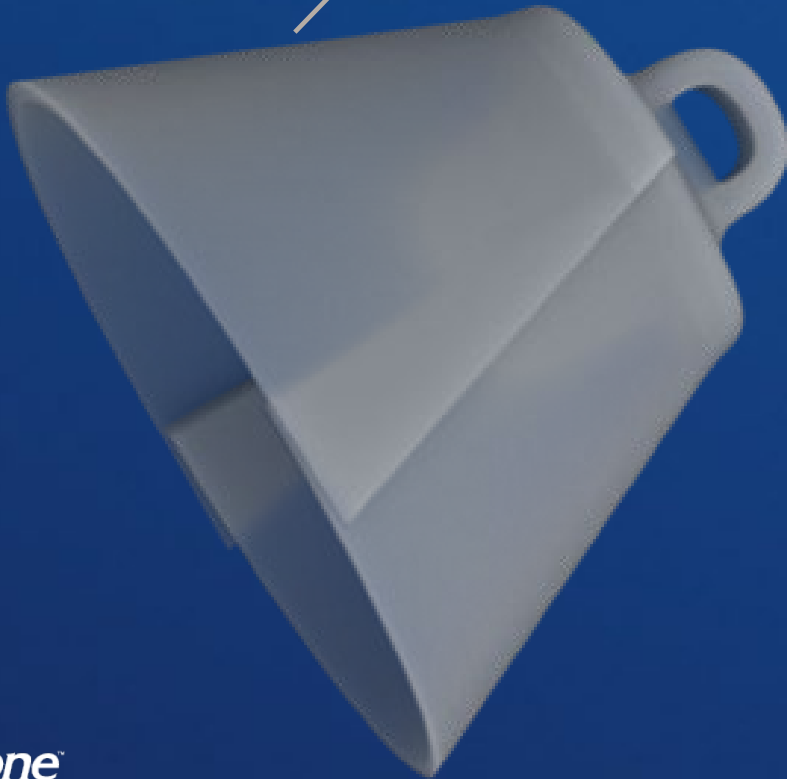

**Tulip dome
design**

2 sizes

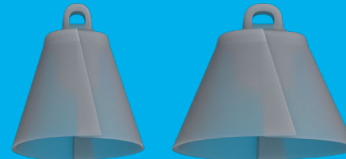


S

M


Tulip Dome

Good choice for a variety of ear canal shapes.

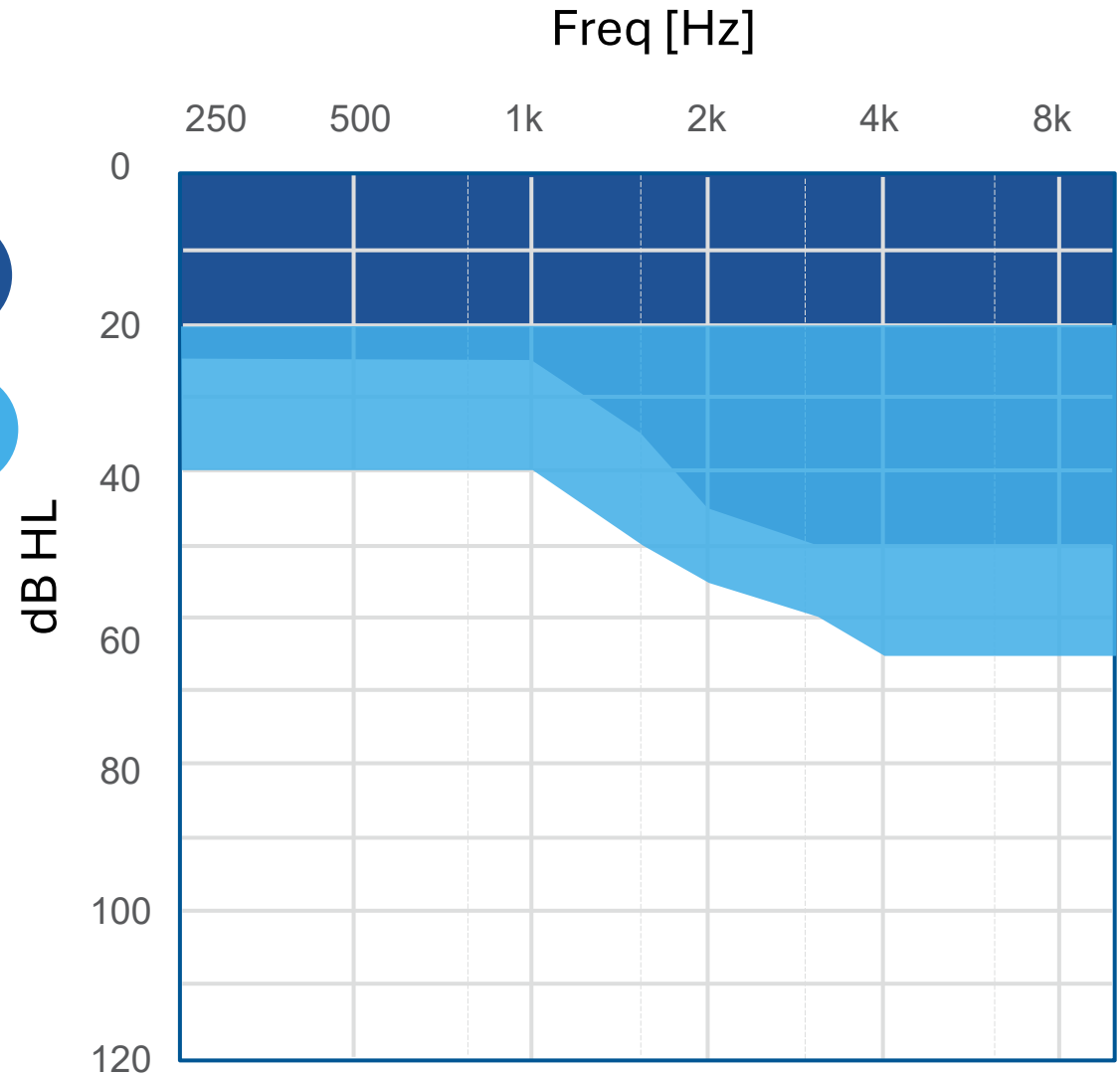

Designed for wearer comfort

Intended to be used as an instant custom-tip style dome and good for narrow ear canals

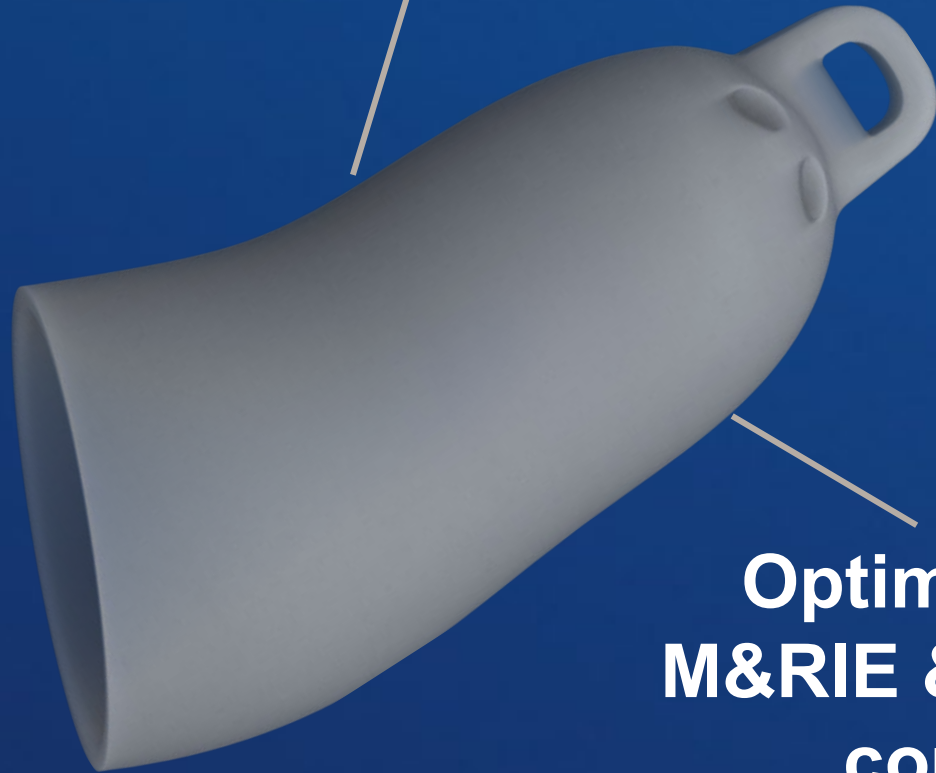

Optimized for
M&RIE & listening
comfort



New Sleeve dome

3 sizes



Sleeve Dome

Provides acoustic comfort of an open dome while providing some occlusion for best use with M&RIE.

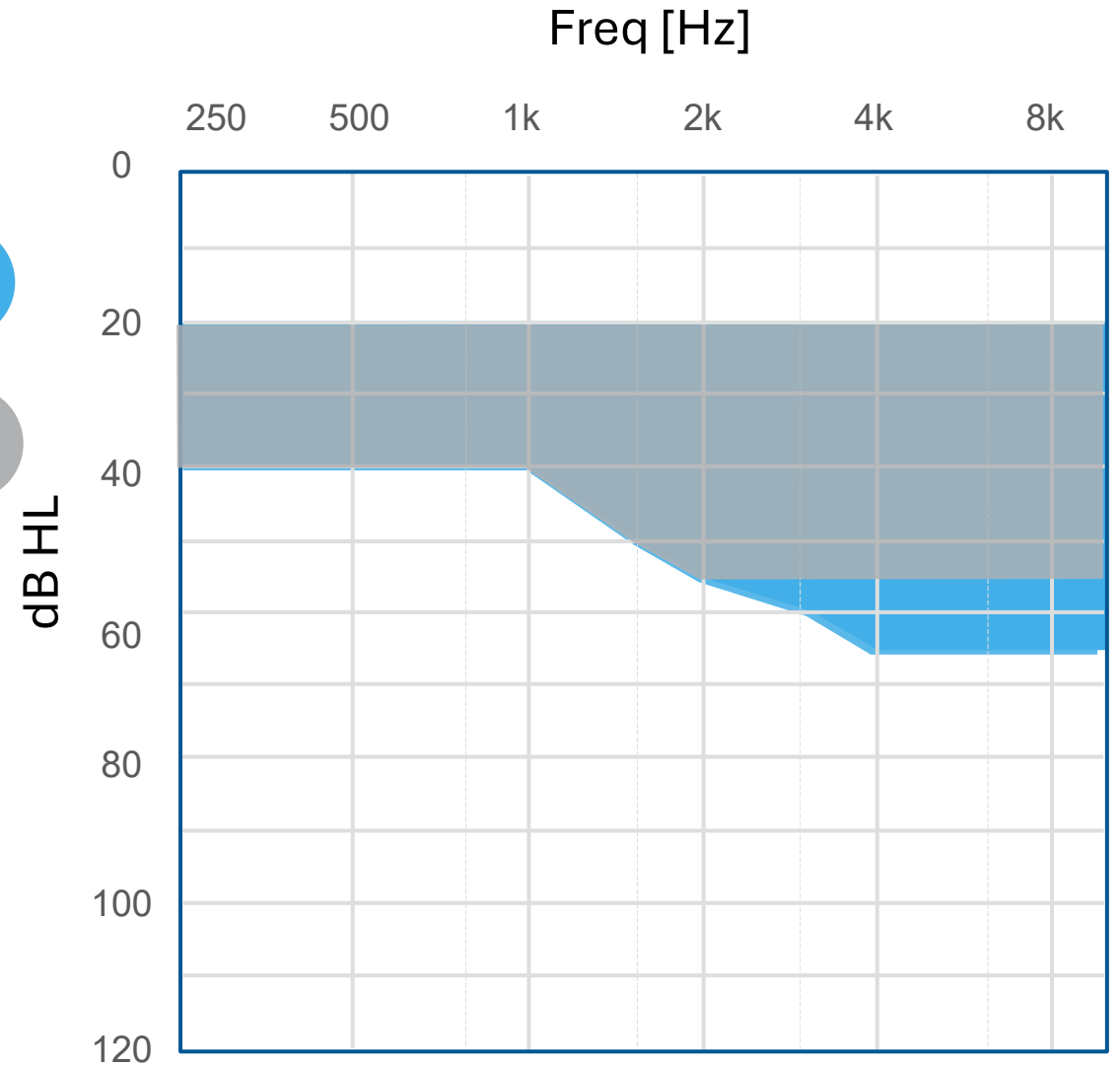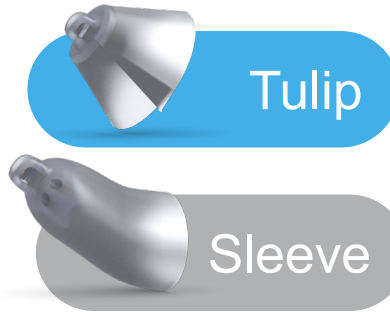

**Additional vent
over legacy**



**New
Vented dome
design**



5 sizes



XS S M L XL

Vented Dome

Offers more occluding fit while still allowing for comfort with its small vents.

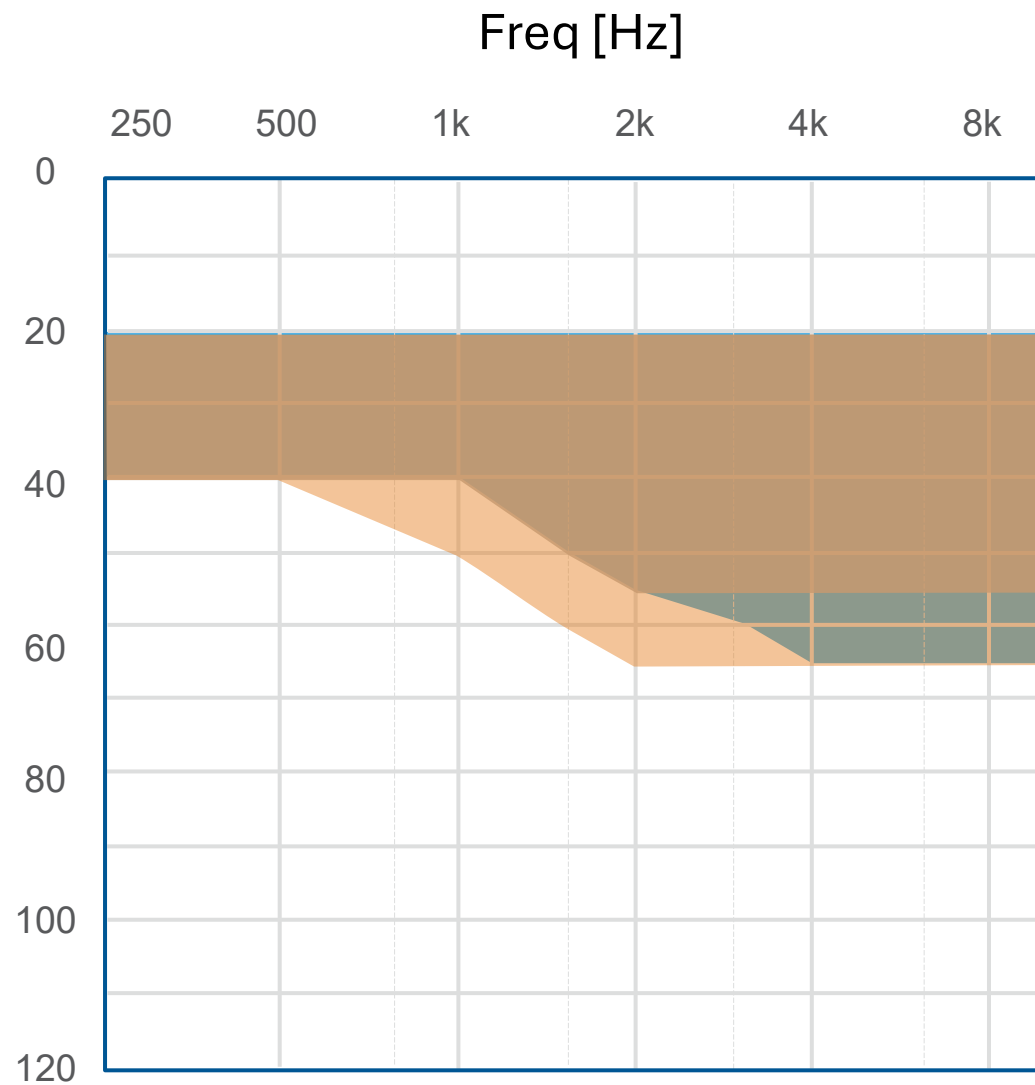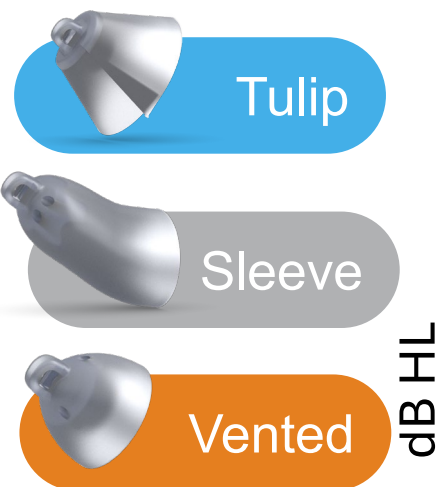

**Superior
streaming
sound quality**



**Power dome
design**

4 sizes



S


M

L


XL

Power Dome


Most occluding dome,
maintains low frequencies
with minimal roll off.

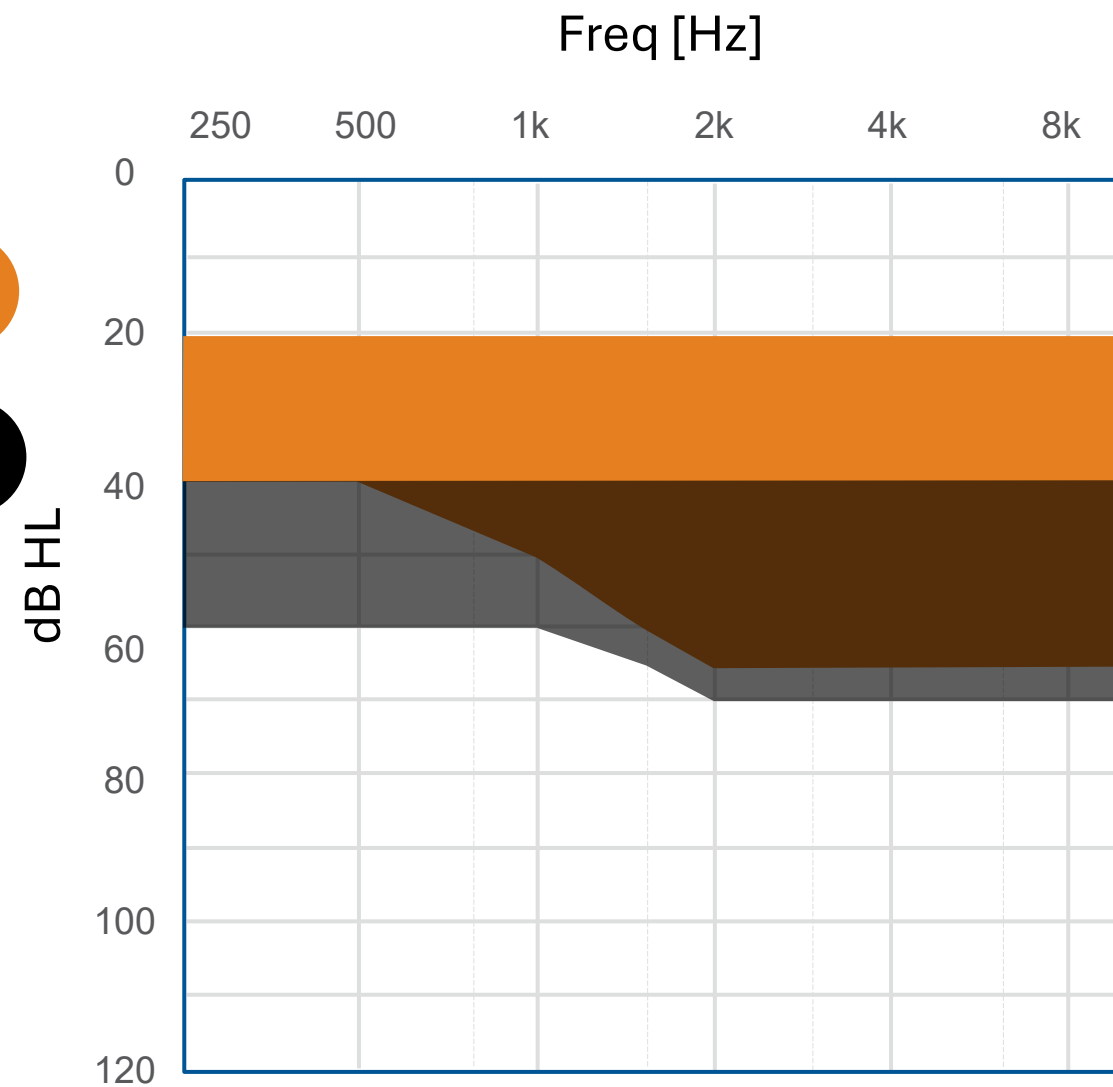

AllFit Dome Frequency Response

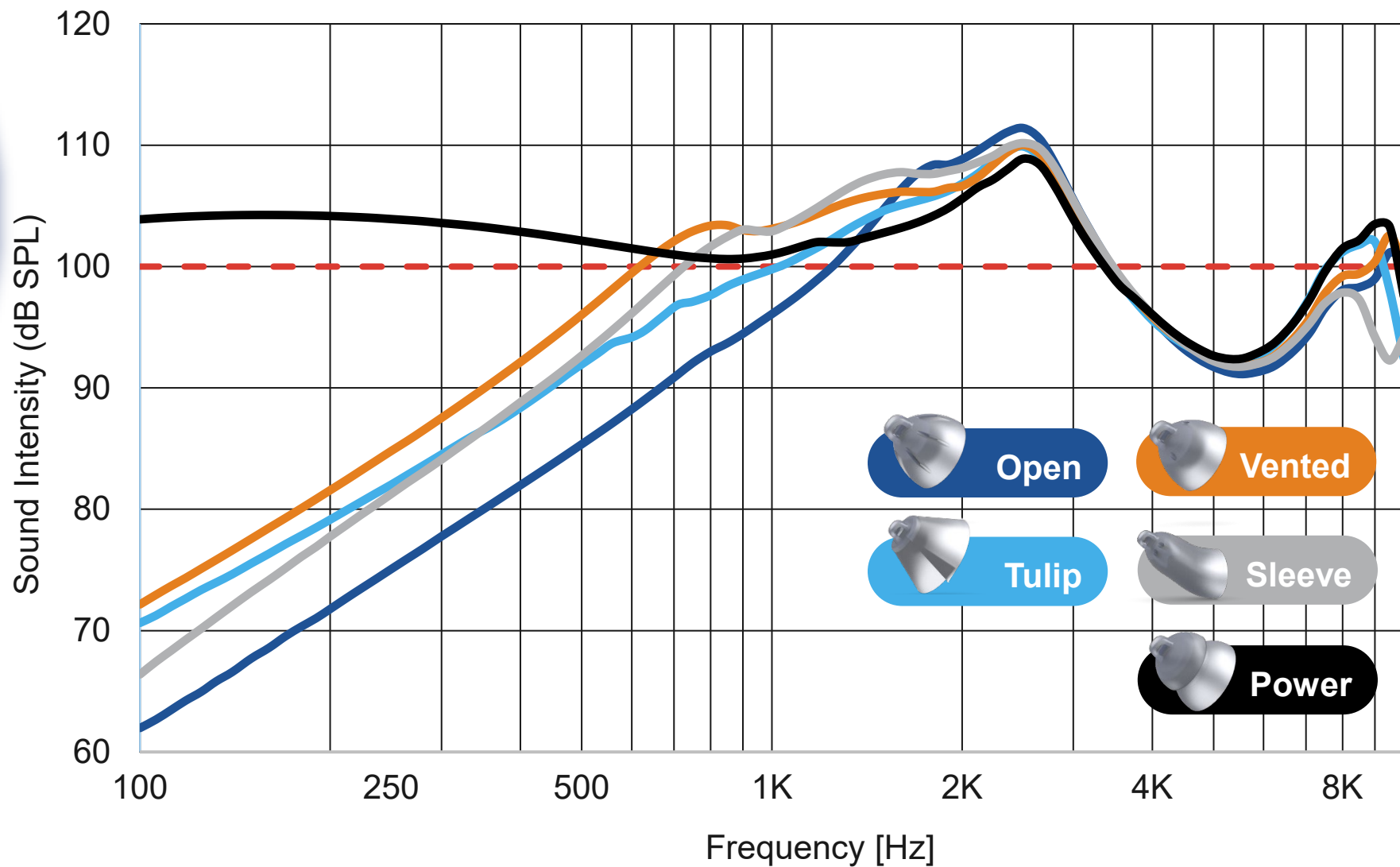

AIIFit Dome Vent Effect

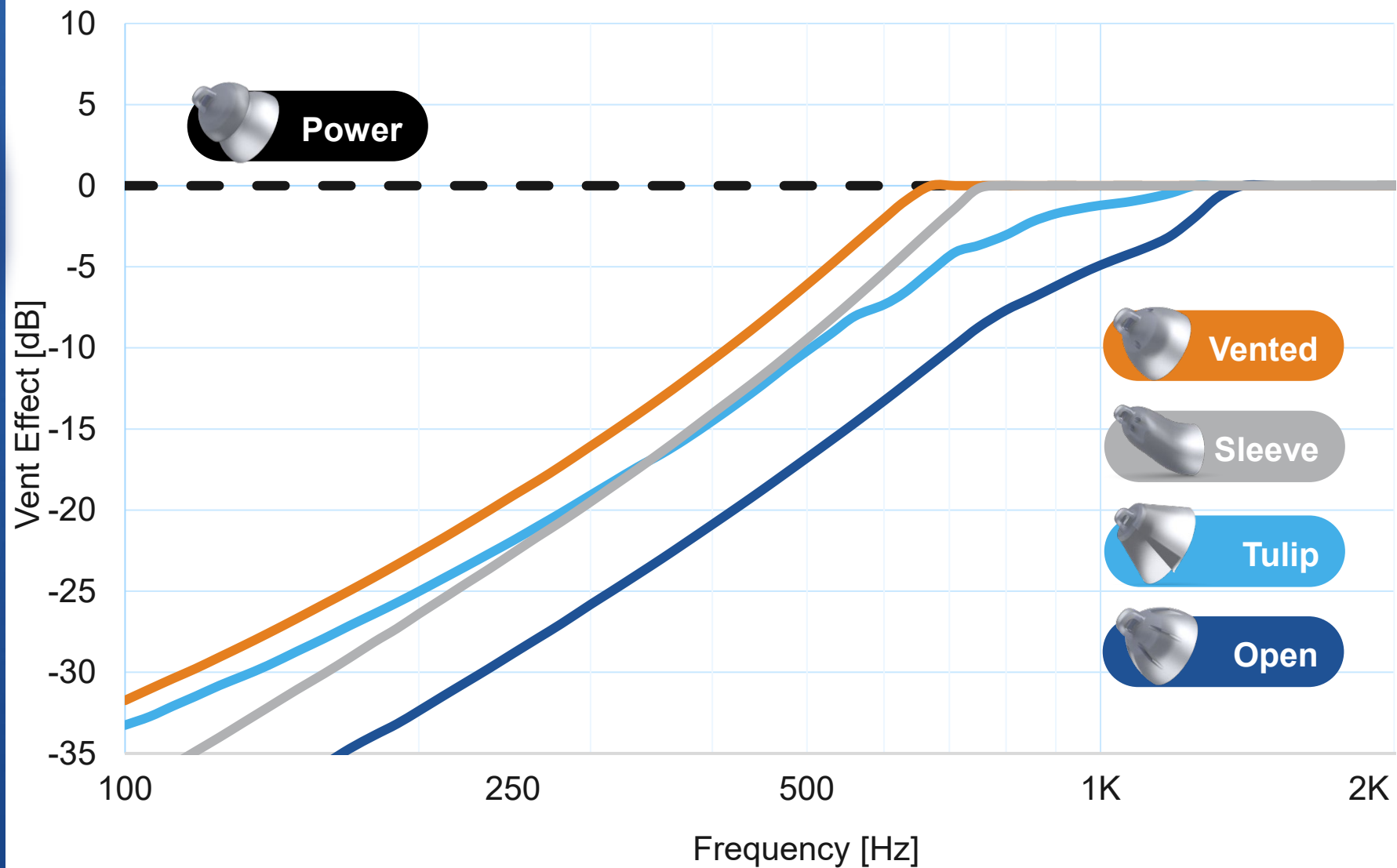

The Correct Fit Provides



Feedback
Prevention



Sound
Delivery



Hearing in
Noise



Audio
Streaming

Wax Protection Loop



NEW

AllFit Domes



More options to fit all ears
More sizes and precise acoustic fits

**Designed to fit
all ears**

97%

Coverage based on
16,000 real ear
impressions¹

**Easier dome
replacement**




**Improved base
response**

Additional wax protection

New tip built into design to
protect better against wax



**Improved
durability**

Increased
silicone thickness



**Designed for best comfort
and acoustic performance**

19 domes



✦ More
comfortable
fit in ear ✦



**Enhanced low
frequency
sound quality¹**

**Sleeve dome
provides better
stable gain, less
opportunities for
feedback**

New sleeve dome
Custom-like fit,
optimized for M&RIE



Learning Outcomes

1

Describe the improvements seen with our new AllFit dome portfolio.

2

List the different occlusion properties for the different AllFit dome styles.

3

Identify the appropriate AllFit dome based on ear canal anatomy, hearing loss configuration, and occlusion considerations.